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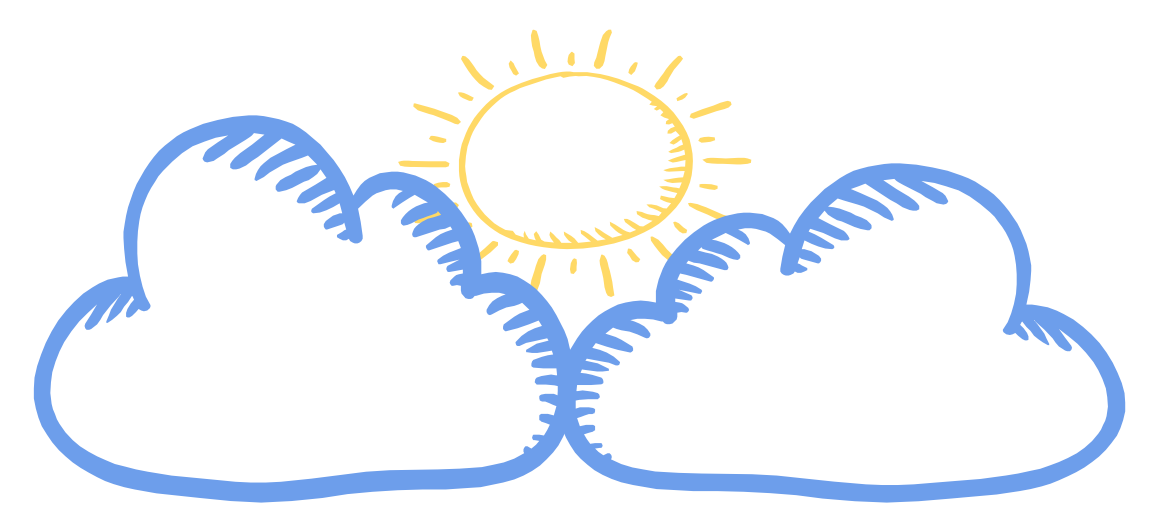
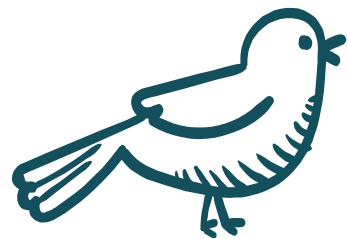
Cochrane Database of Systematic Reviews

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
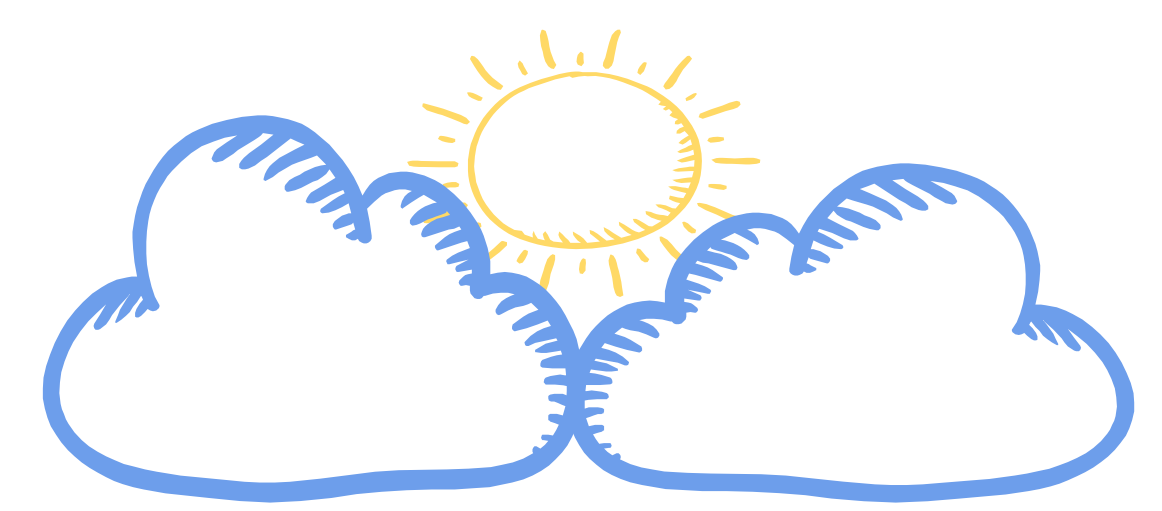
Cochrane Database of Systematic Review 2013, Issue 7.

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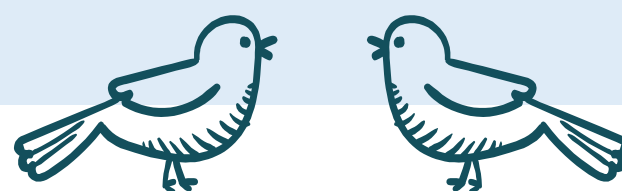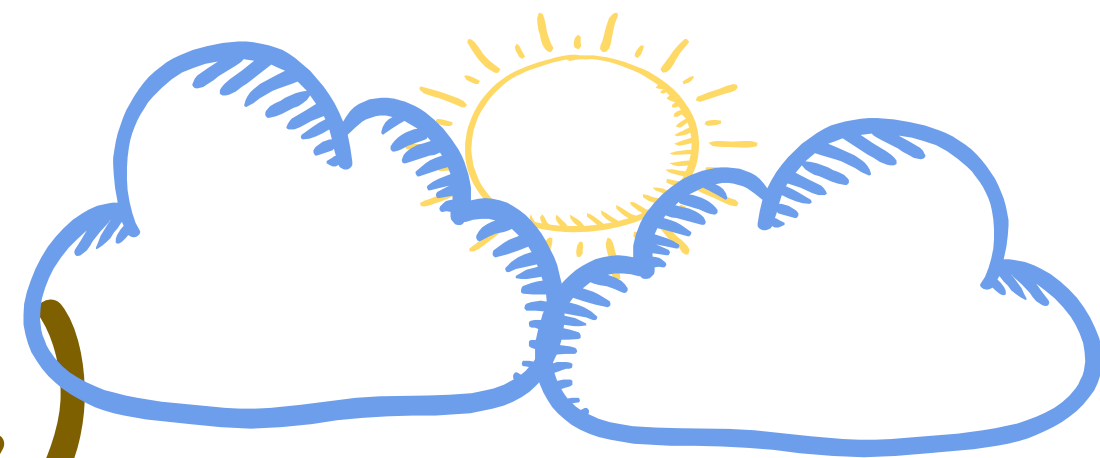
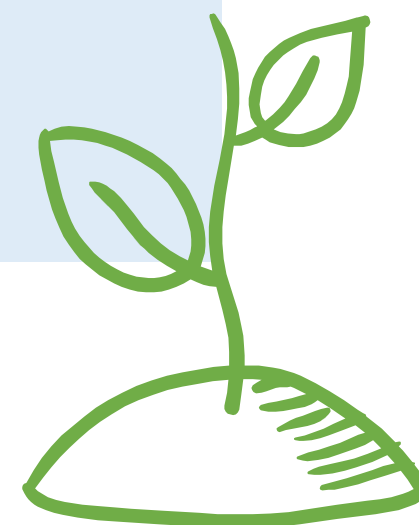
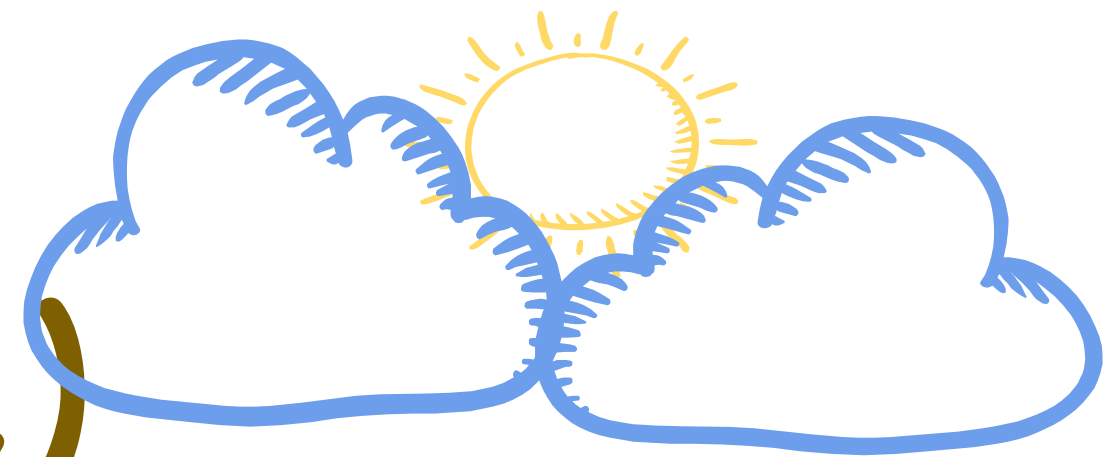
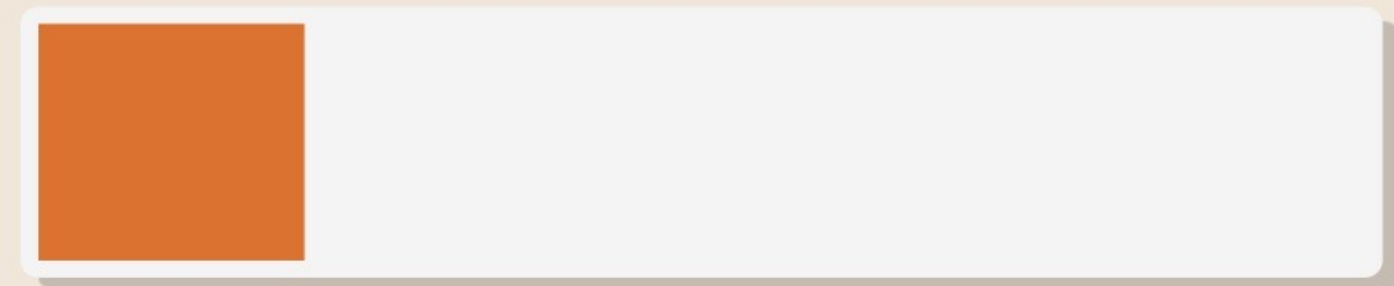
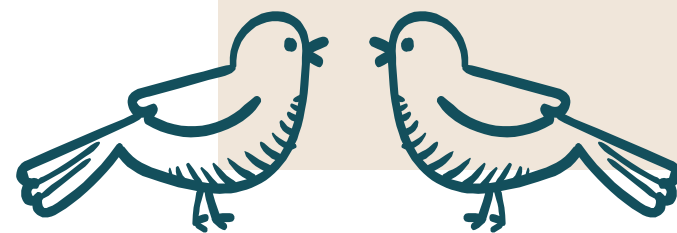
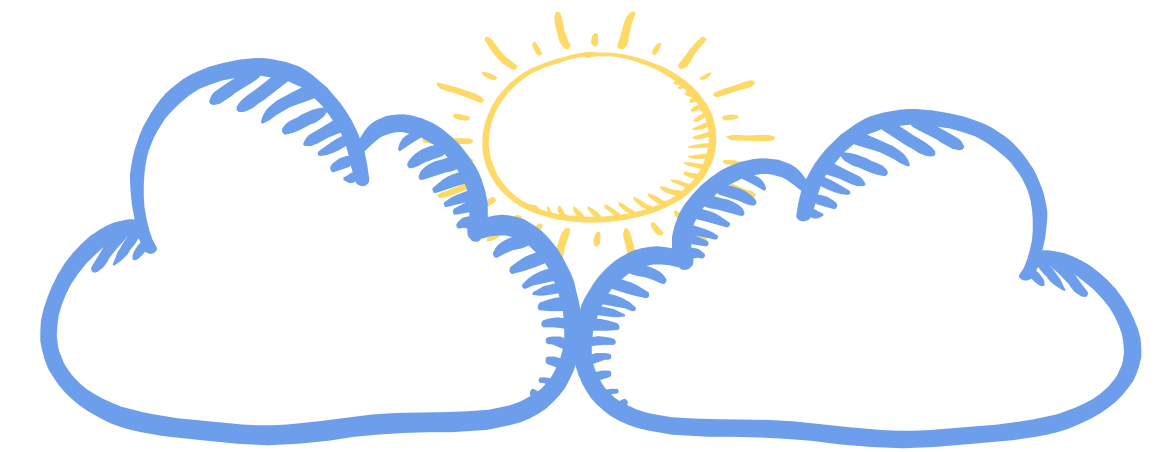
	Fluoride varnish
Fluoride ion	22,600 ppm.
Amount of use	0.3-0.5 ml.
Swallowing of fluoride	6.8-11.3 mg
Remain of fluoride in Enamel	longer
Advantage and Disadvantage	can apply for specific tooth surface



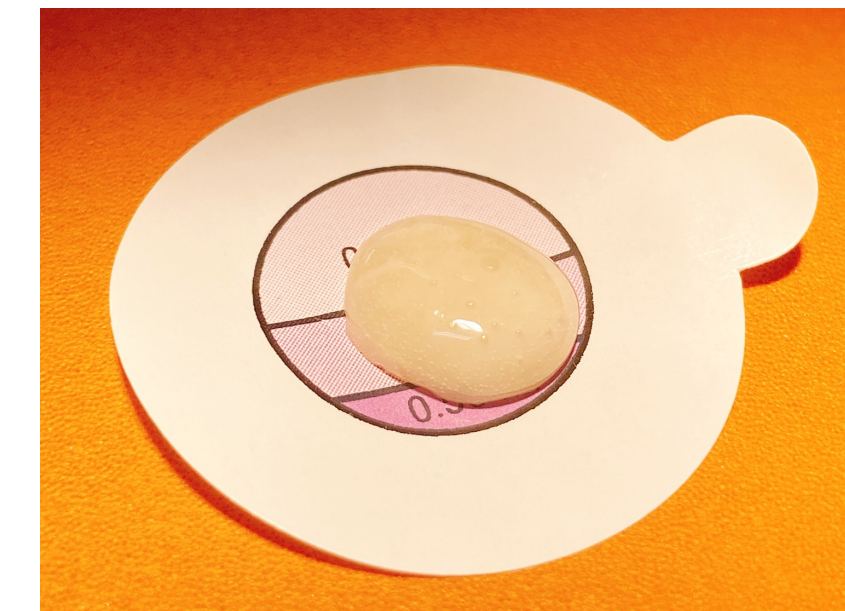
Photo by Dr.Worawit Sakulthai

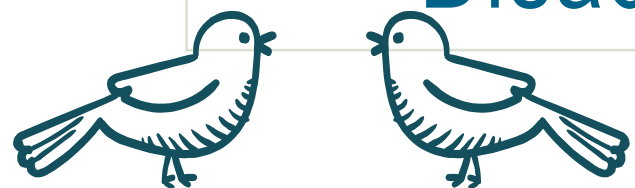Fluoride Gels


Fluoride gels

- Consist of 1.23% acidulated phosphate fluoride (APF) or 2% neutral sodium fluoride (NaF)
- Need to be worn in a tray for four minutes for maximum effect
- Not recommended for patients ages six and under due to a higher risk of swallowing the fluoride.¹
- Time-consuming to set up for dental professionals
- Uncomfortable for patients²

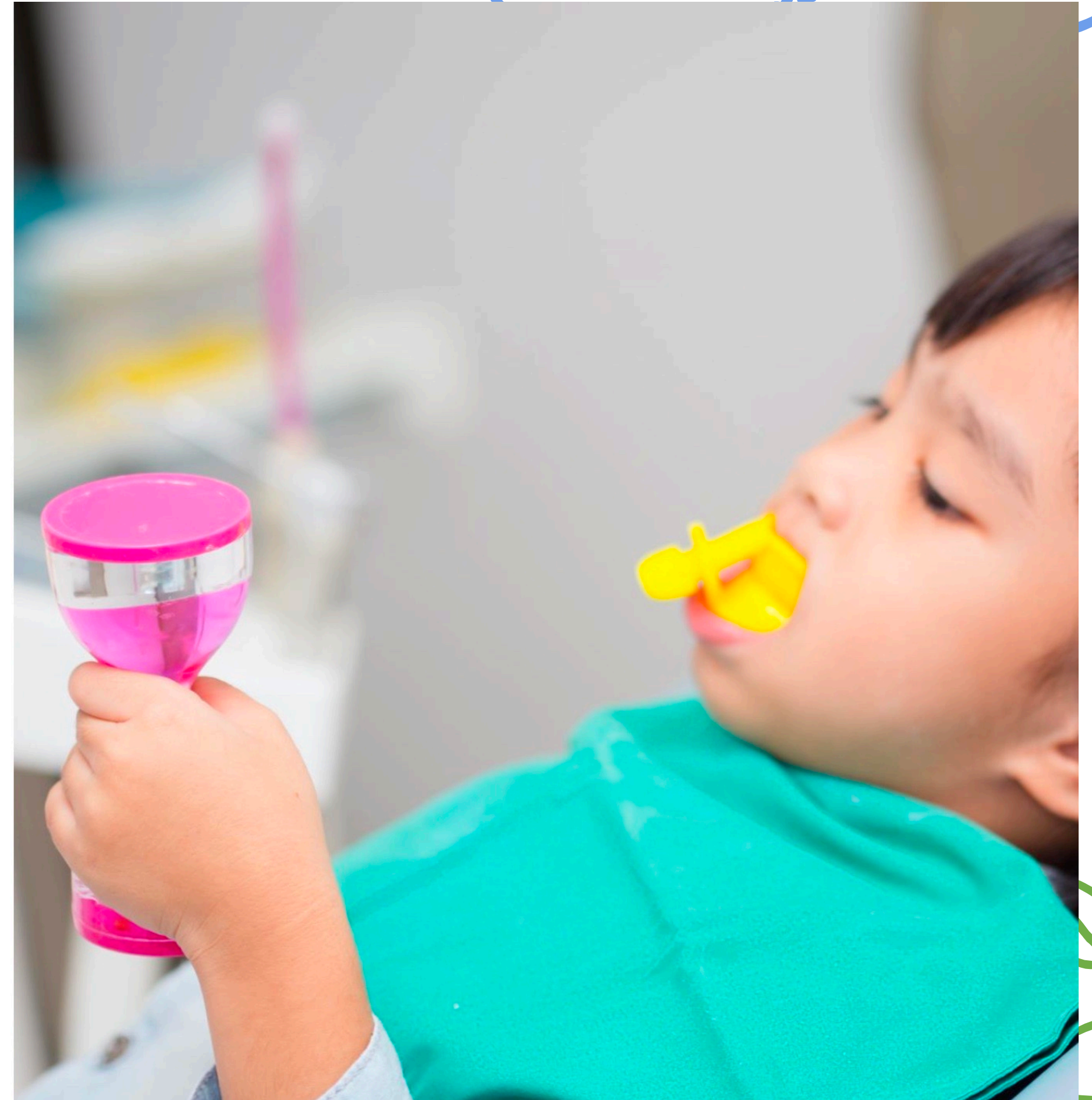
1. JADA 2013;144(11):1279-1291
2. R. Hawkins et al. J Public Health Dent 2004;64(2): 106-110

Fluoride varnish vs. fluoride gel

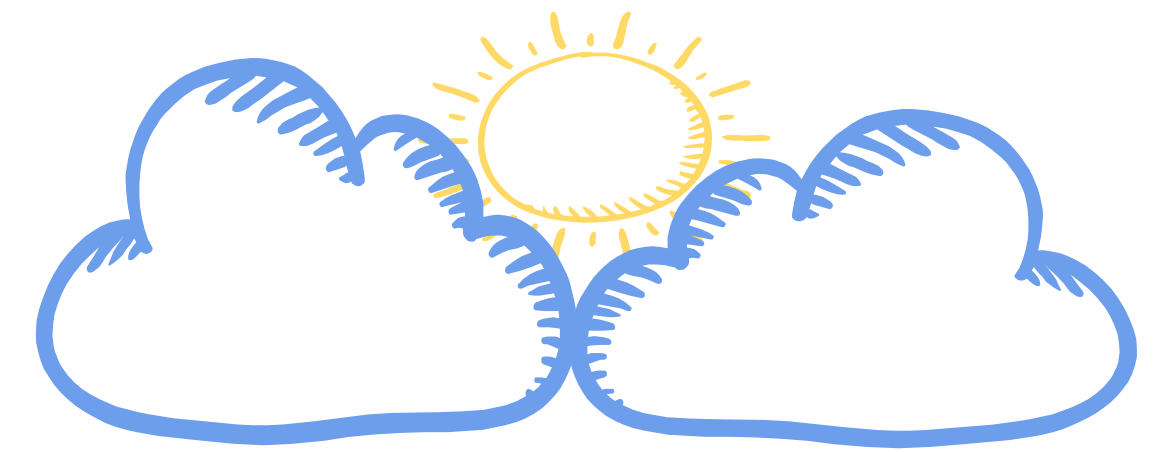	Fluoride varnish	Fluoride gel
Fluoride ion	22,600 ppm.	12,300 ppm.
Amount of use	0.3-0.5 ml.	5 ml.
Swallowing of fluoride	6.8-11.3 mg	10-35 mg (without suction) 2-7 mg (with suction)
Remain of fluoride in Enamel	longer	shorter
Advantage and Disadvantage	can apply for specific tooth surface	do not use in child <3 years old

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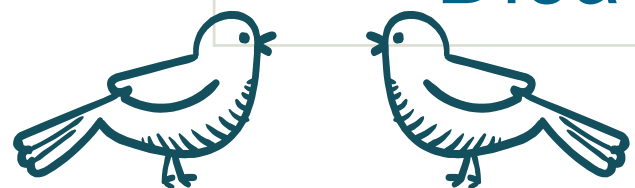Less-than-ideal experiences

Traditional fluoride treatments deliver unpleasant experiences for both patients and professionals

Foams, gels, and rinses¹

- Multi-step process
- Uncomfortable tray
- Higher risk of swallowing
- No hygienic unit dose
- Unpleasant taste


Fluoride varnish

- Mixing required
- Sticky
- Unpleasant taste, smell, and mouthfeel²


1. R. Hawkins et al. Prevention. British Dental Journal 2003; 195: 313–317
2. Source: Internal. Clinical Study Report.

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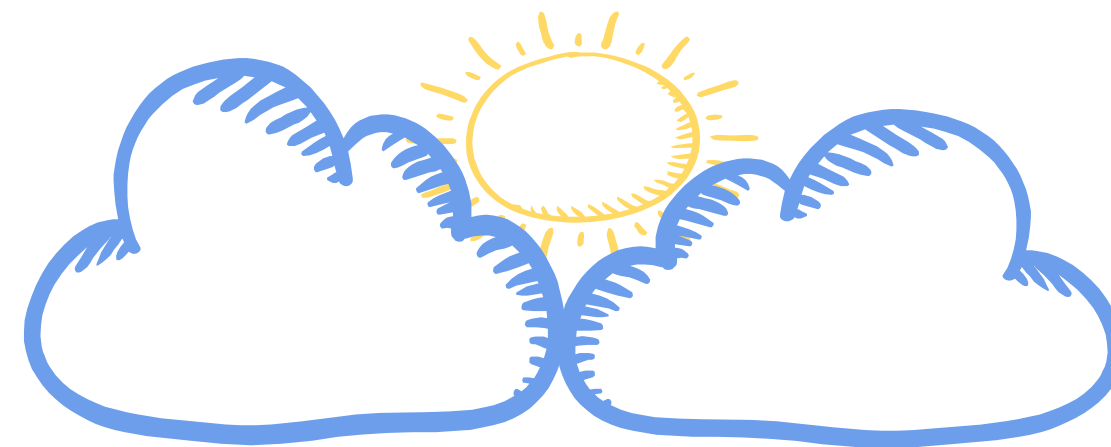
Fluoride	% caries reduction	
	permanent	primary
Water fluoridation	50–60	40–50
Fluoridated milk	35–77	40–69
Fluoride tablet	20–40	50–58
Fluoride mouth rinse	23–30	
Fluoride toothpaste	23–32	
Fluoride gel	18–25	
Fluoride varnish	30–46	7.4–43.8

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




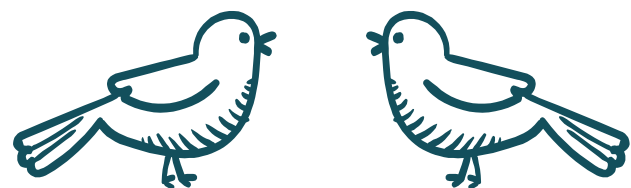
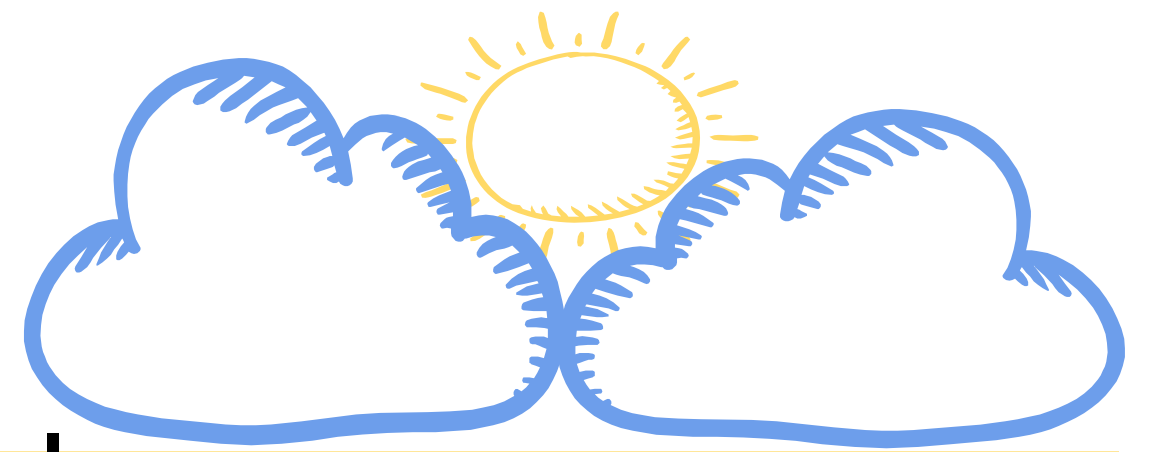
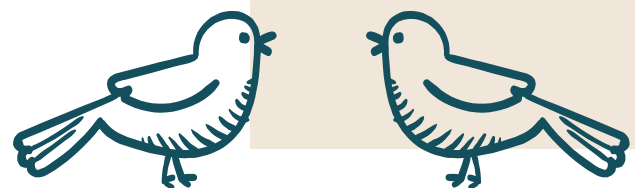
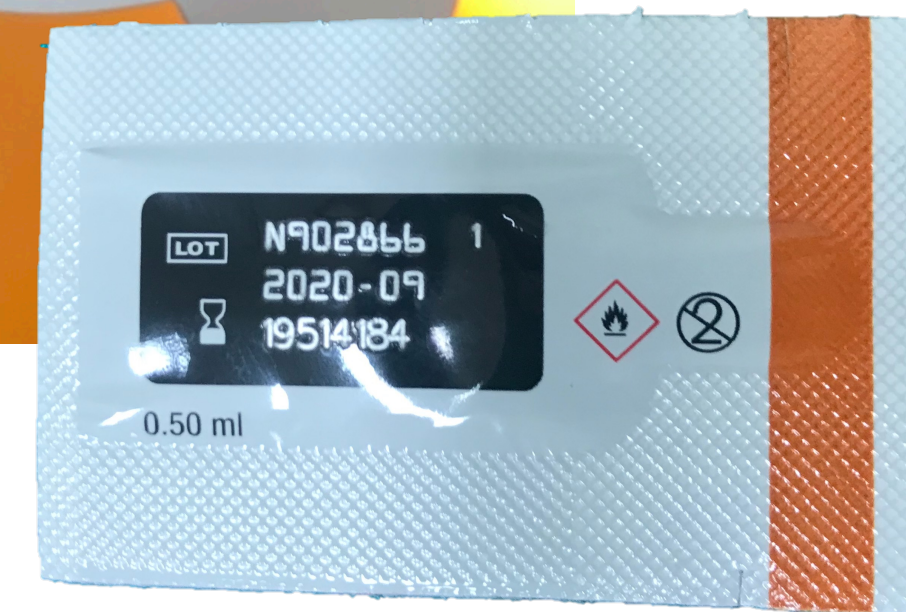
Photo by Dr.Worawit Sakulthai

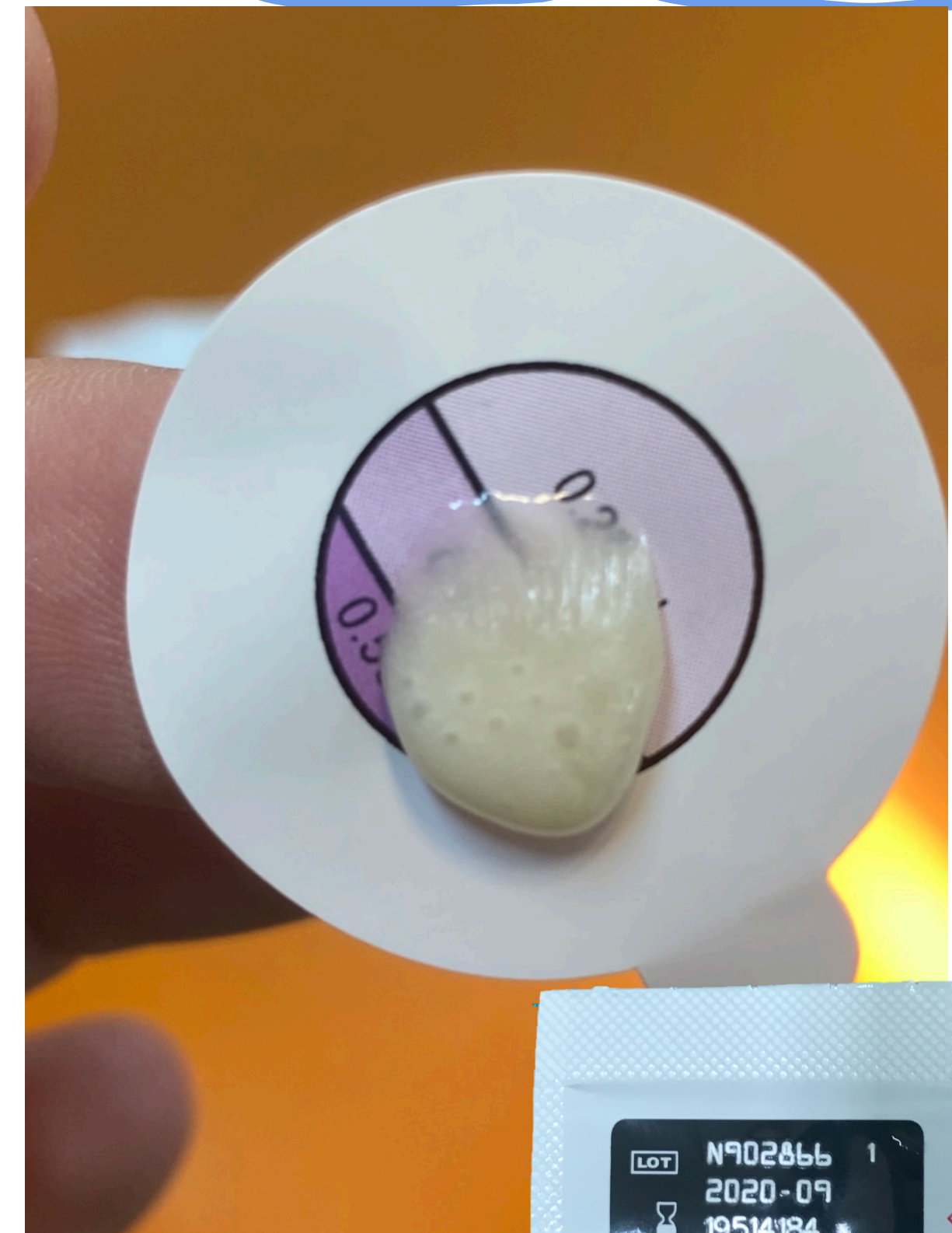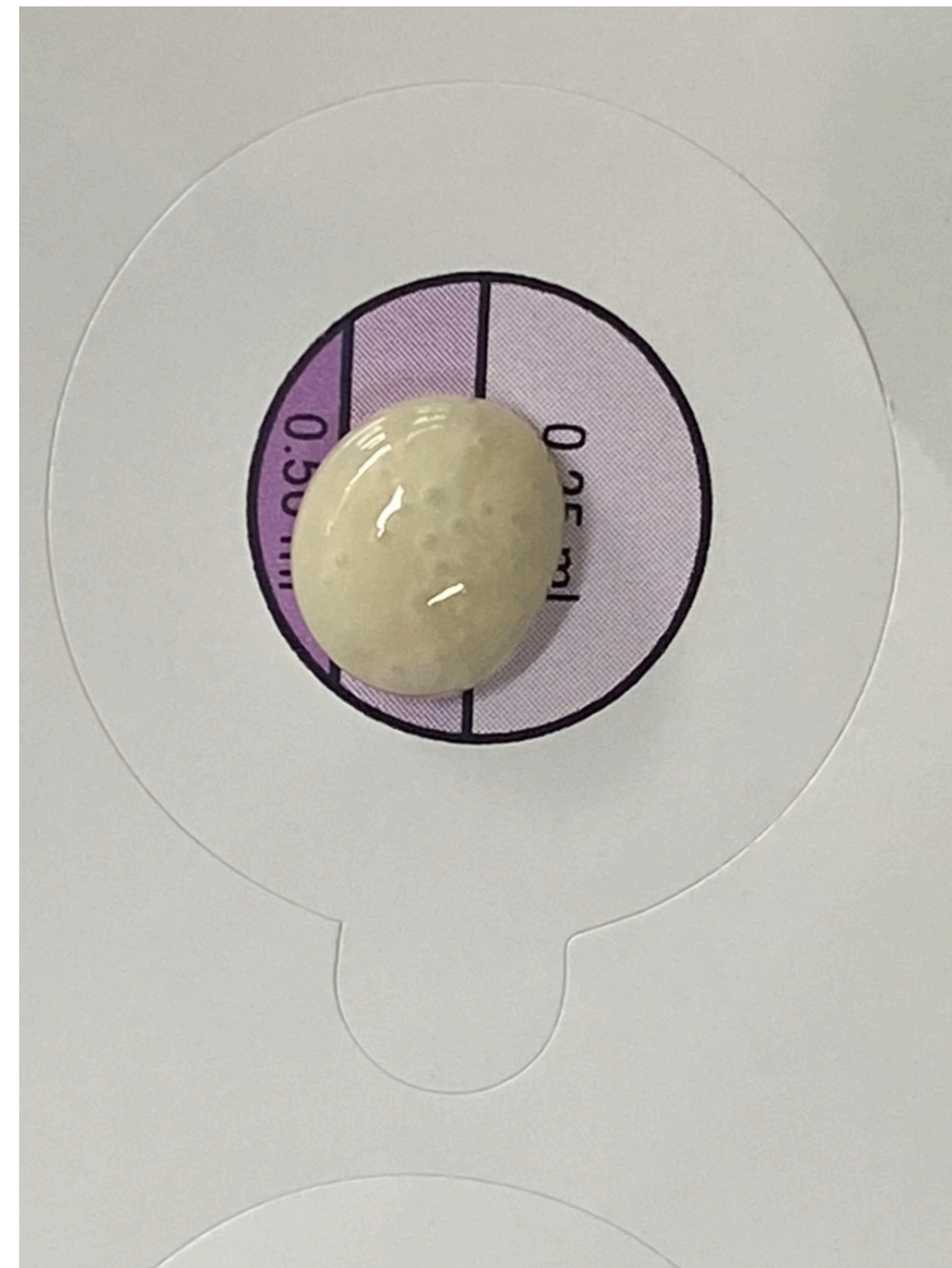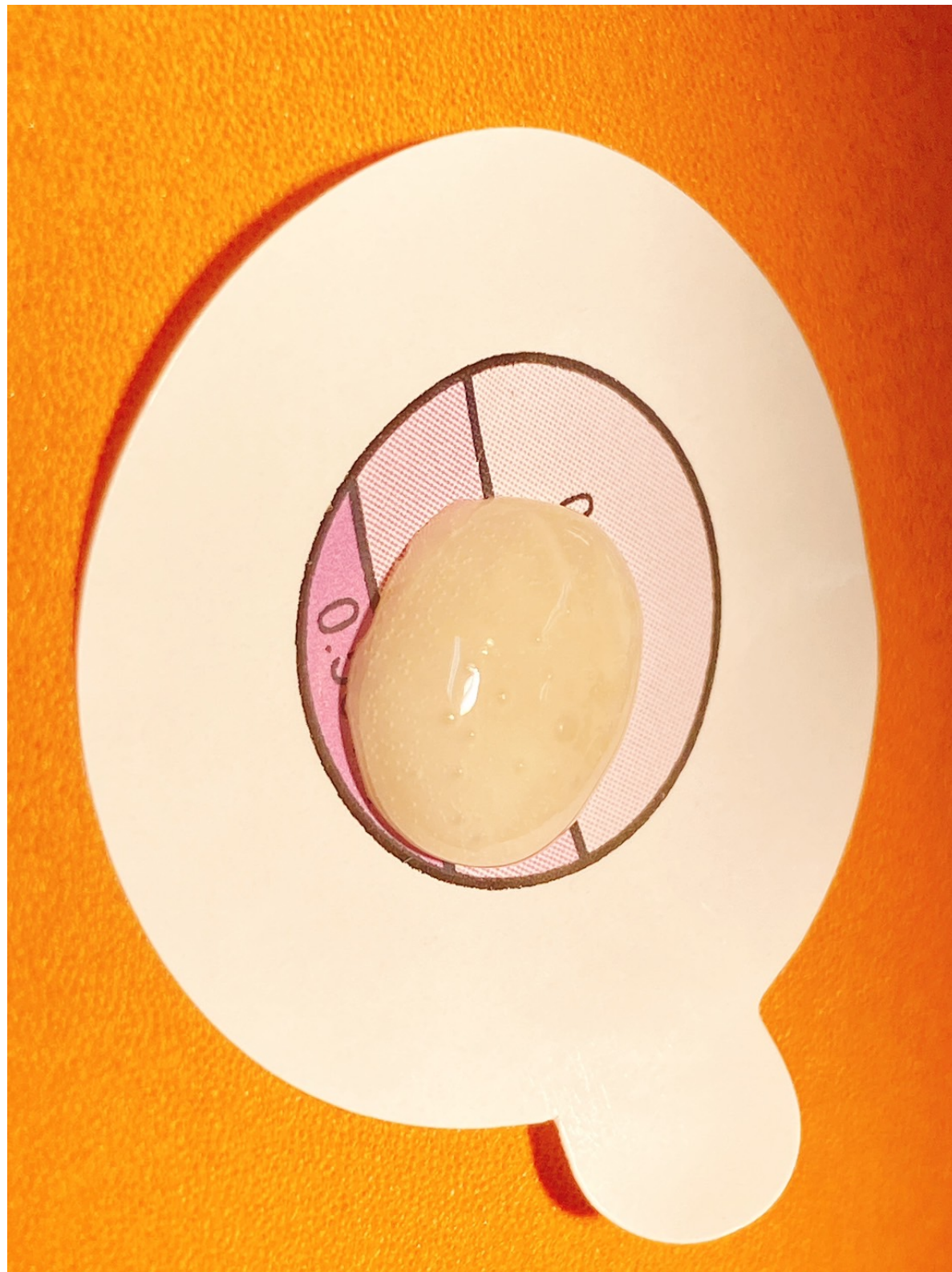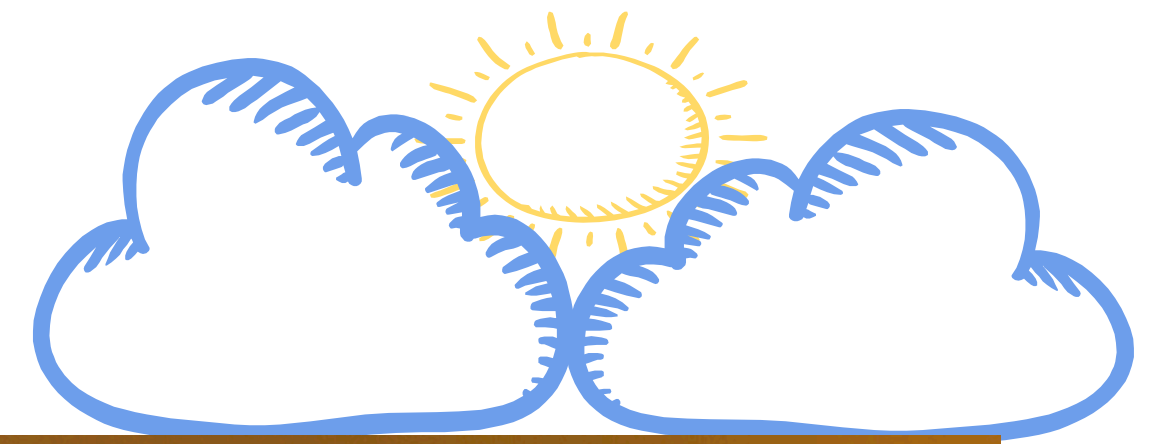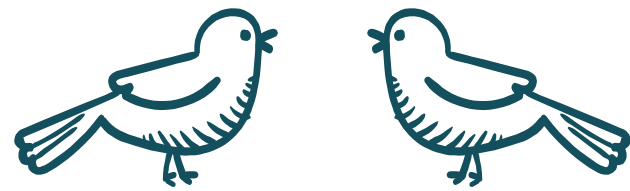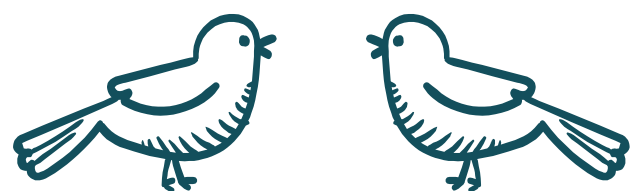Fluoride varnish

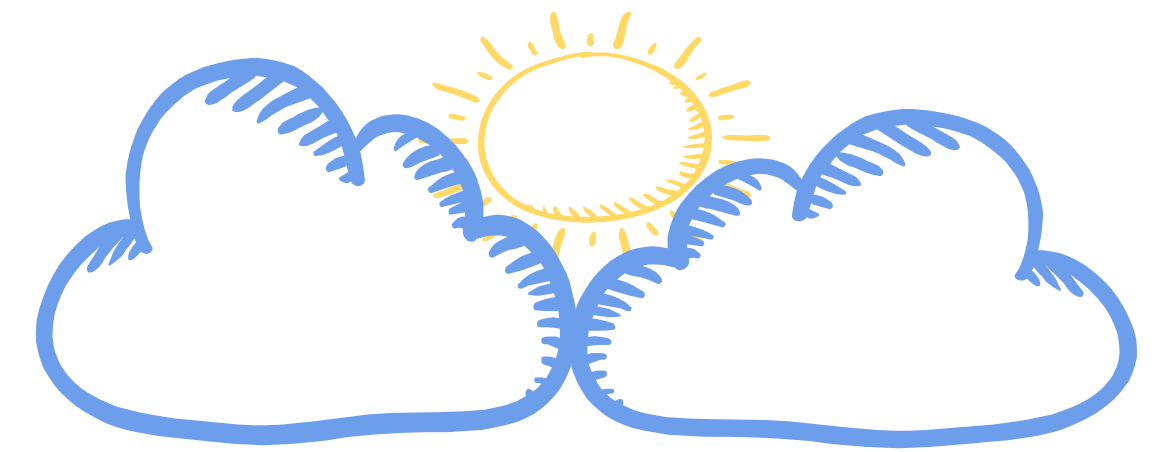
Traditional fluoride varnish treatments

- Consist of rosin, solvents, and 5% sodium fluoride (NaF) in suspension
- Varnish must be stirred right before application
- Much like oil and water, they don't mix
- Rosin traps fluoride in place, requiring a longer treatment time for the ions to be dissolved
- Varnish is sticky, adhering to teeth and other surfaces
- When a fluoride varnish contacts saliva, it creates a clumpy, tacky, rough texture on teeth
- Uncomfortable feeling can last up to 24 hours

Photo by Dr.Worawit Sakulthai


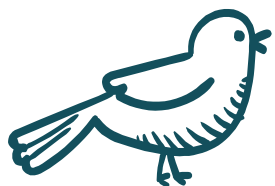
Fluoride varnish

Fluoride's Mechanism of Action

1. Enhancing Remineralization

- Fluoride promotes the uptake of calcium and phosphate from saliva to demineralized enamel.
- Forms **fluorapatite** ($\text{Ca}_{10}(\text{PO}_4)_6\text{F}_2$), which is **more resistant to acid dissolution** than hydroxyapatite.

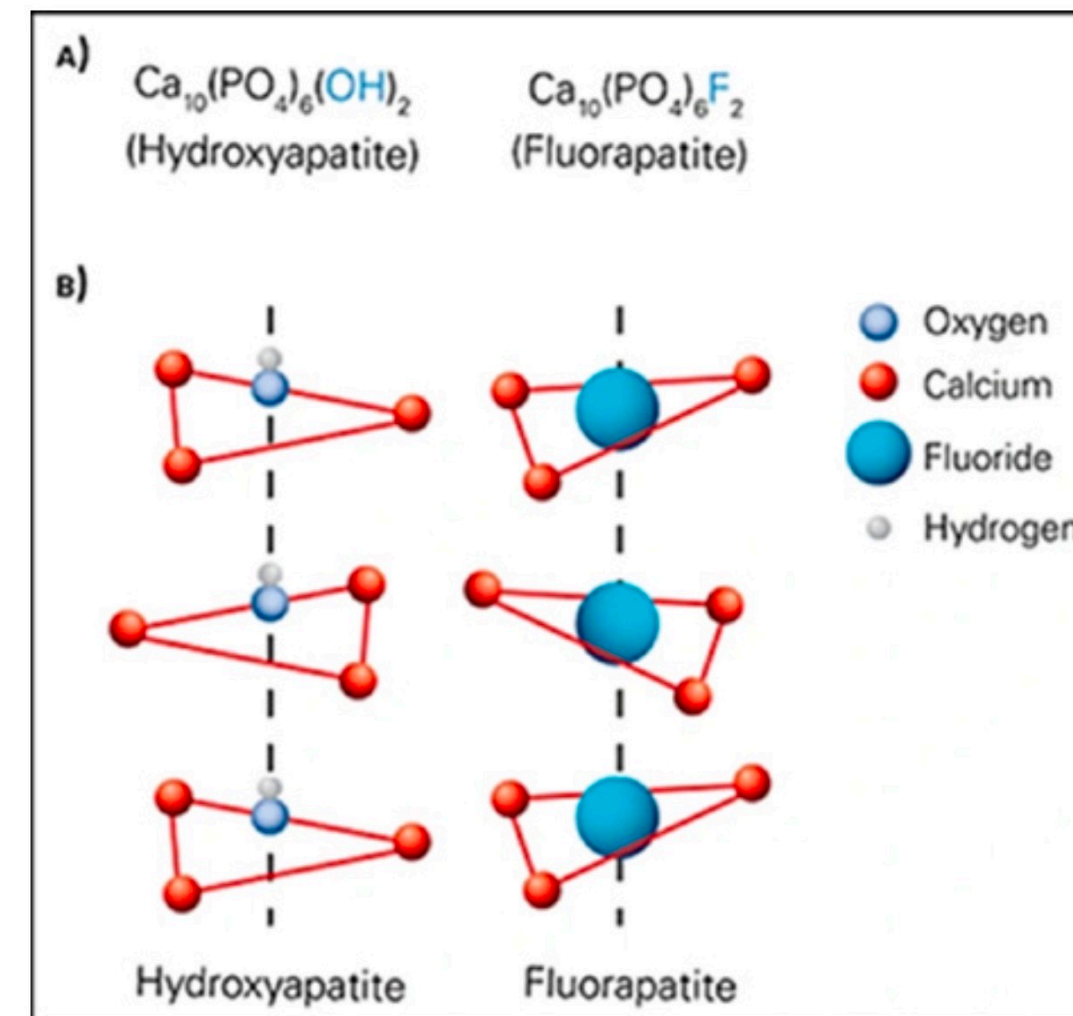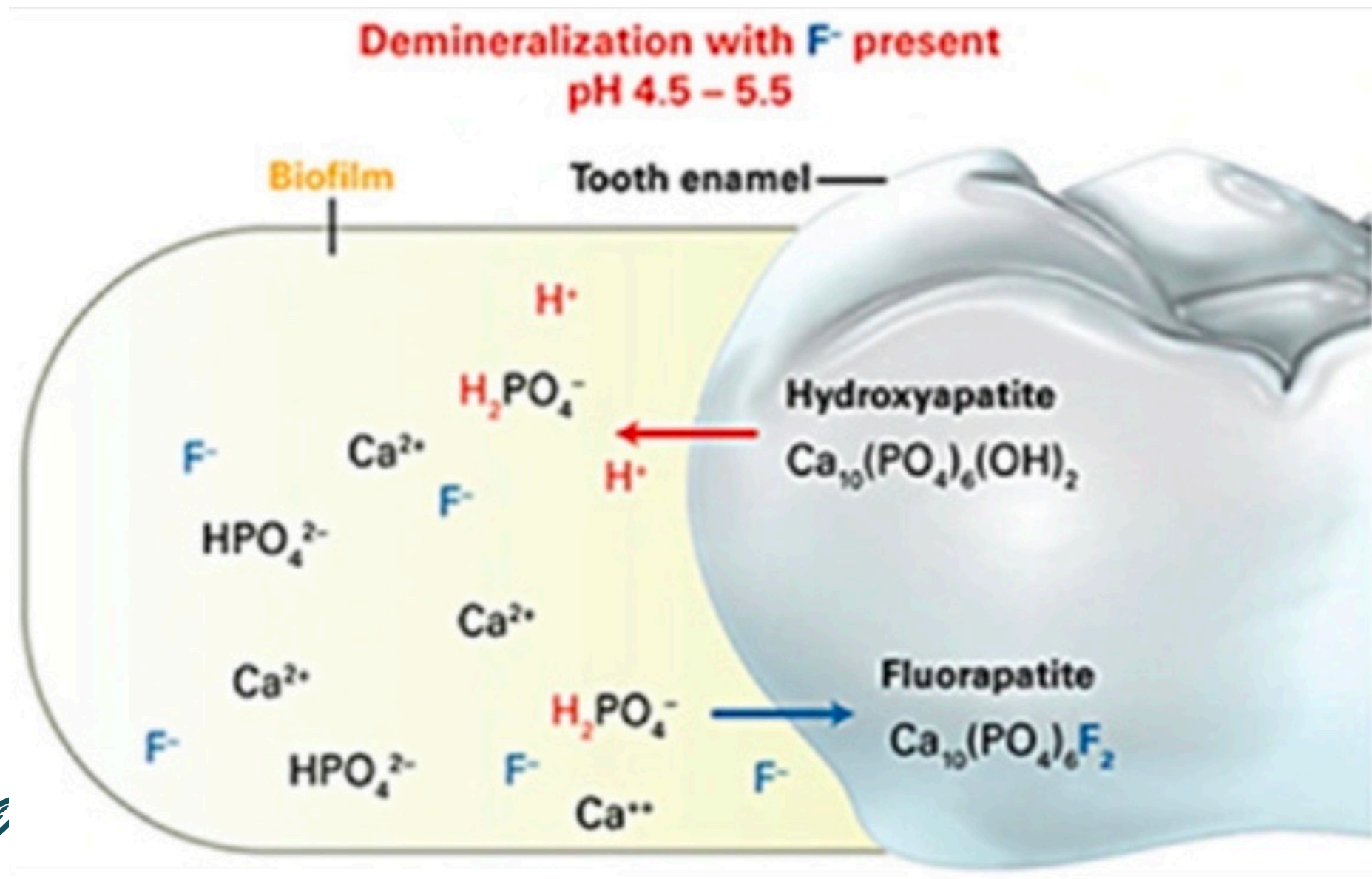
Fluoride varnish


Fluoride's Mechanism of Action

2. Inhibiting Demineralization

- When fluoride is present in plaque or saliva, it **adsorbs to the enamel surface** and forms a protective barrier.
- This lowers the solubility of enamel during acid attacks.

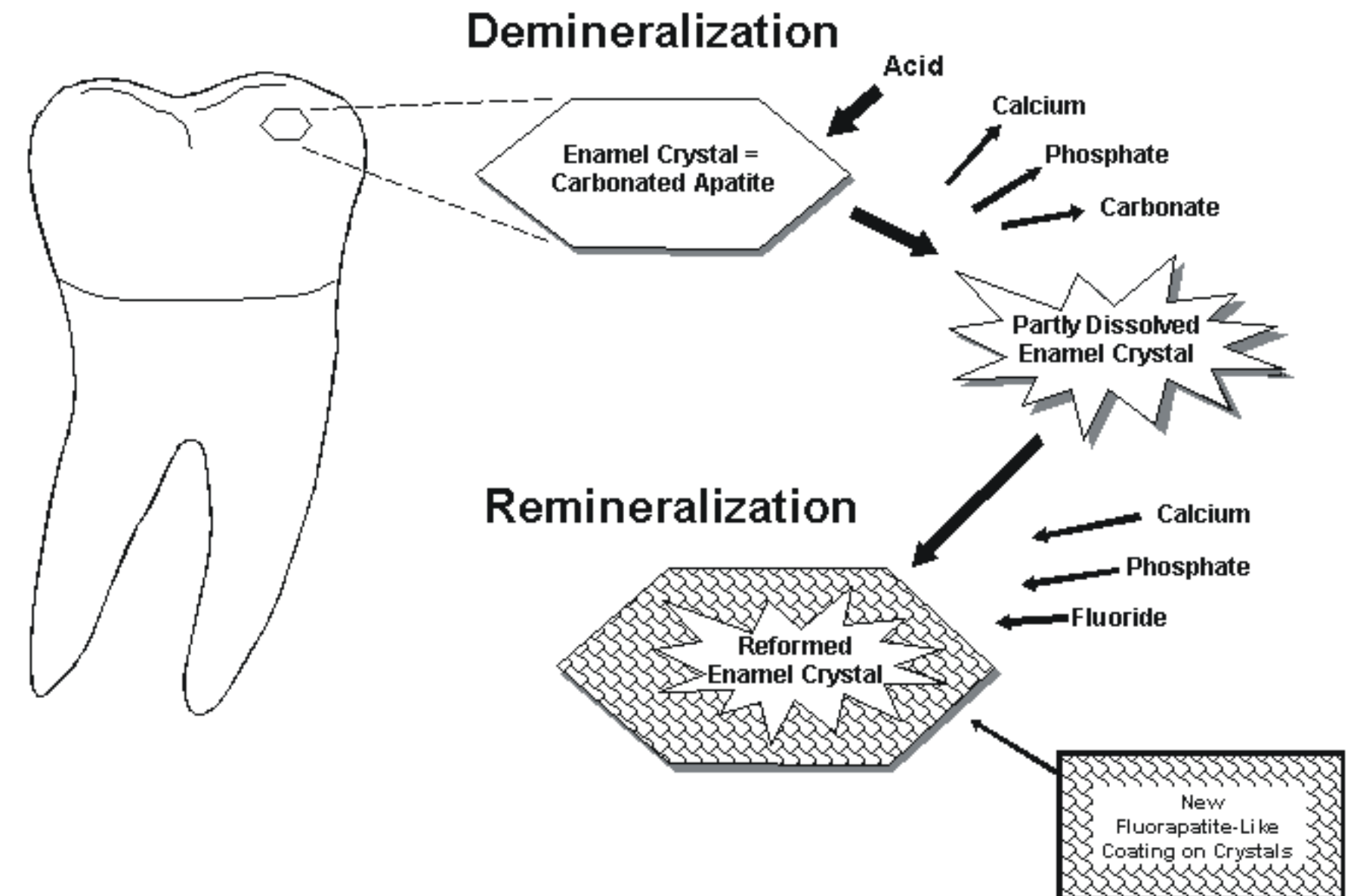William Maas et al, 2019

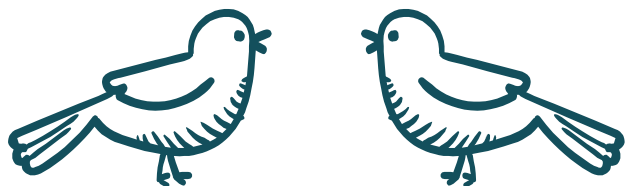
Fluoride varnish

Fluoride's Mechanism of Action

3. Antibacterial Effects

- At high concentrations (e.g., fluoride varnish), fluoride is **bactericidal**, especially to *S. mutans*.
- At low concentrations (e.g., water fluoridation), it is **bacteriostatic**—inhibiting key enzymes like enolase, reducing acid production.

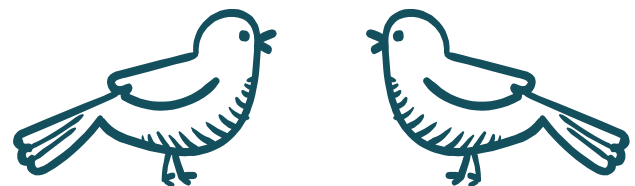
Fluoride varnish


Fluoride's Mechanism of Action

4. Formation of CaF_2 -like Reservoirs

- Topical fluoride application forms **calcium fluoride-like globules** on the tooth surface.
- These act as **reservoirs**, slowly releasing fluoride during acid challenges.

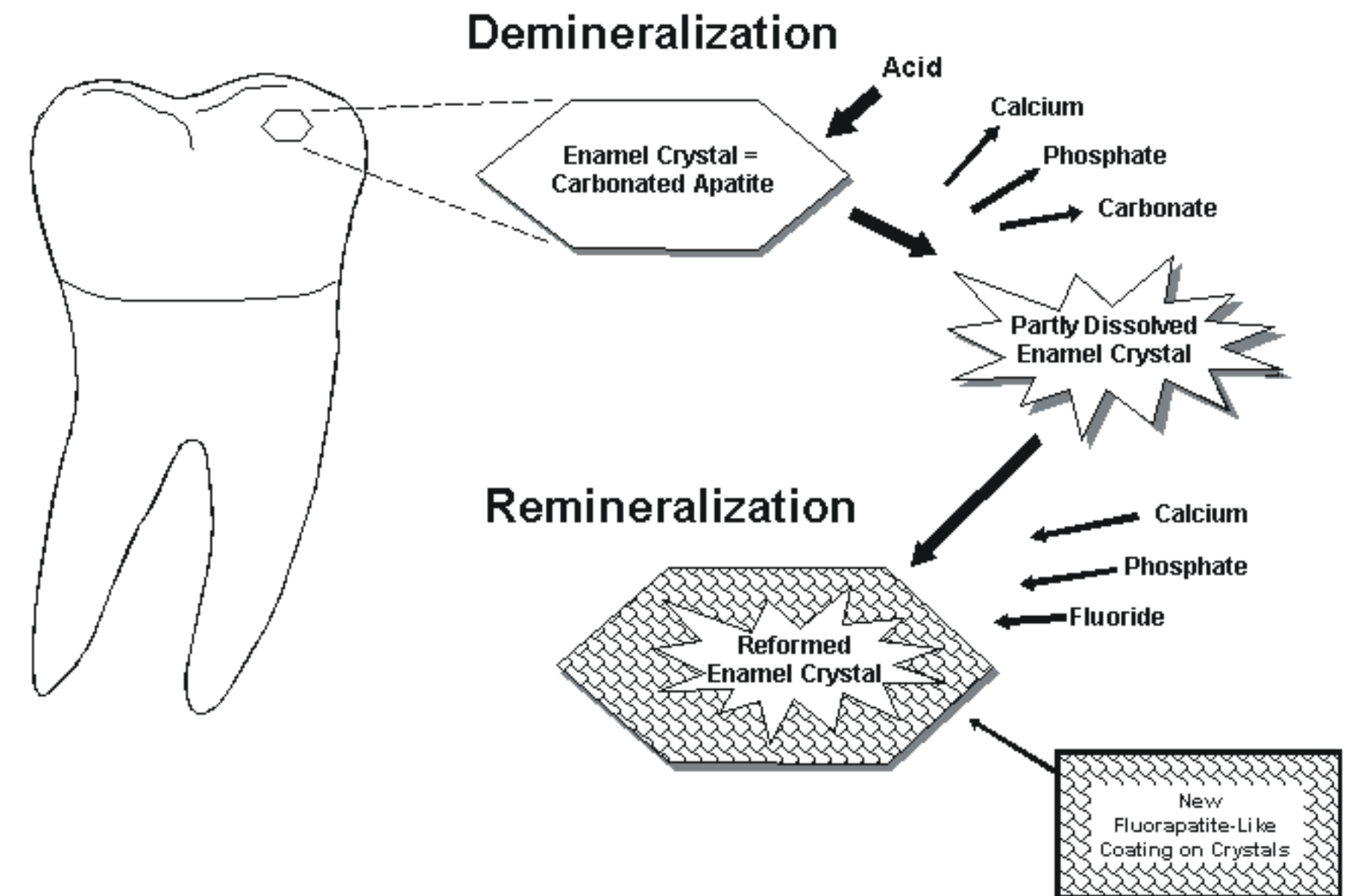William Maas et al, 2019

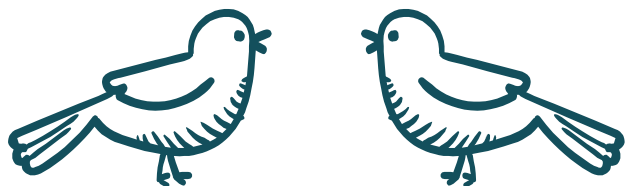Fluoride varnish

3M™ Clinpro™ Clear Fluoride

Water-based technology provides a smoother application

Fluoride varnish

3M™ Clinpro™ Clear Fluoride

Product formulation:

- Water-based formula
- Fluoride ion exists in solution, fully dissolved and ready to release
- Thixotropic (shear thinning) rheology for easy product handling
- Adheres to the tooth surface
- Stabilized calcium phosphate technology
 - Stabilizes fluoride, calcium, and phosphate, preventing premature interaction
 - Enables release of soluble and bioavailable ions
- Contains xylitol
- Flavorless, mint, watermelon

Fluoride varnish

3M™ Clinpro™ Clear Fluoride

Fluoride. Simplified.

3M™ Clinpro™ Clear Fluoride Treatment is the next generation in fluoride treatment

- 2.1% sodium fluoride (NaF)
- Water-based with neutral pH
- Formula free of common allergens*
- 3M stabilized calcium phosphate technology
- Unique, 0.5 ml L-Pop delivery system
- No mixing—just press and go
- Goes on quickly and easily
- Clear appearance on teeth
- 15-minute contact time
- Available in flavorless, mint, and watermelon

*See allergen disclosure at [3M.com/ClinproClear](https://www.3M.com/ClinproClear)

Fluoride varnish

3M™ Clinpro™ Clear Fluoride

No mixing required


Unmixed traditional rosin-based varnishes

- Must be mixed to disperse heterogeneous parts

3M™ Clinpro™ Clear Fluoride Treatment

- Homogenous, water-based solution means no mixing is needed
- Goes on clear and stays that way


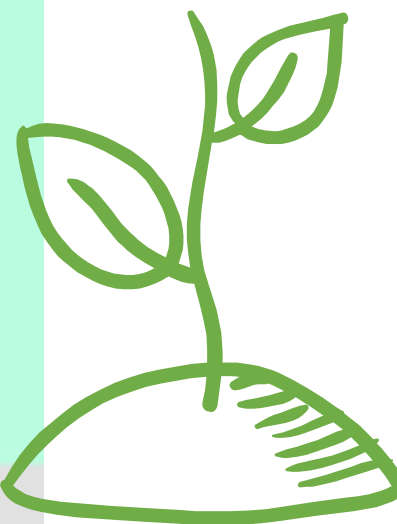
Fluoride varnish

3M™ Clinpro™ Clear Fluoride

Enamel remineralization

Clinpro Clear Fluoride Treatment helps remineralize enamel lesions which can aid in the prevention of caries

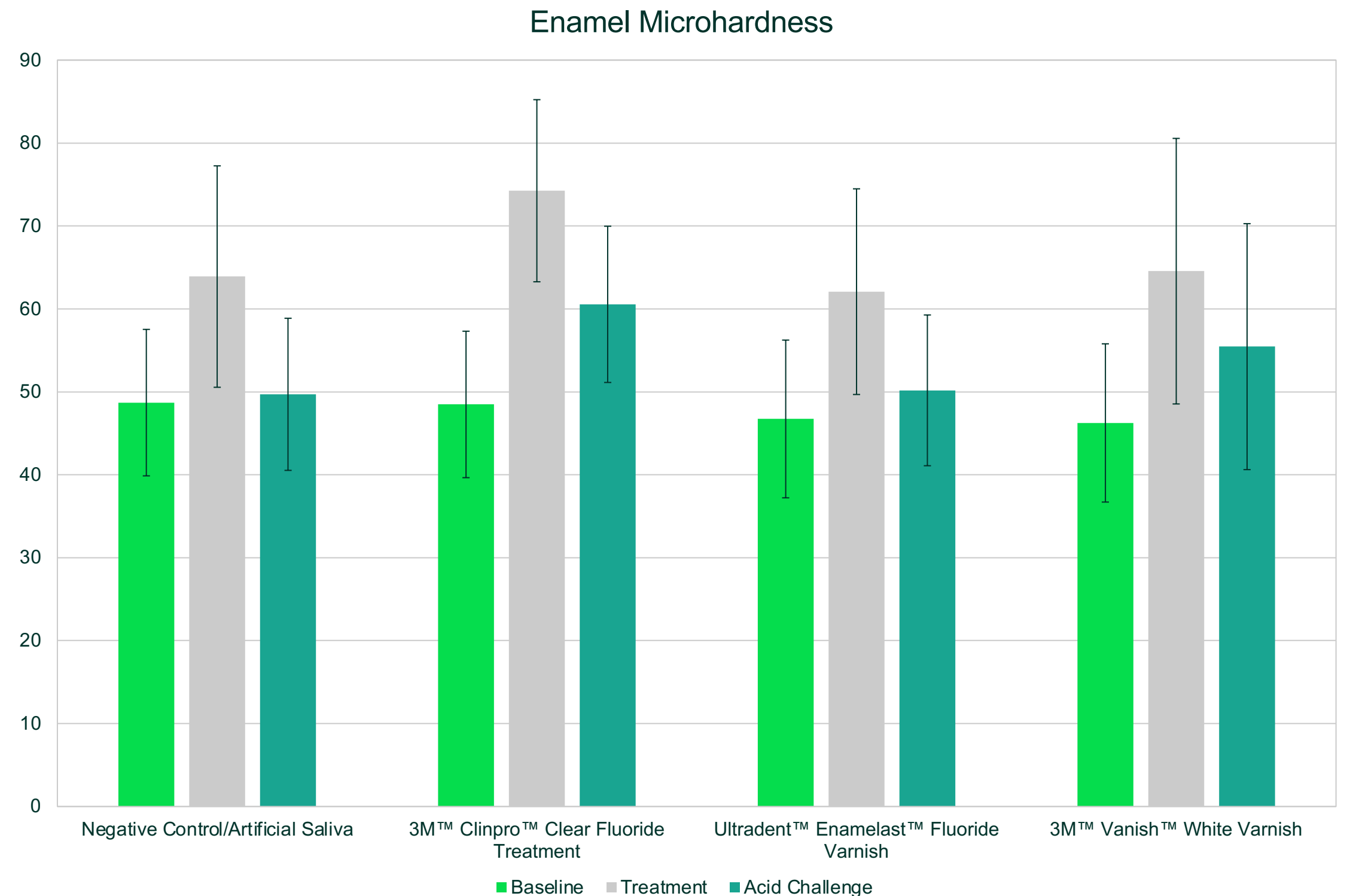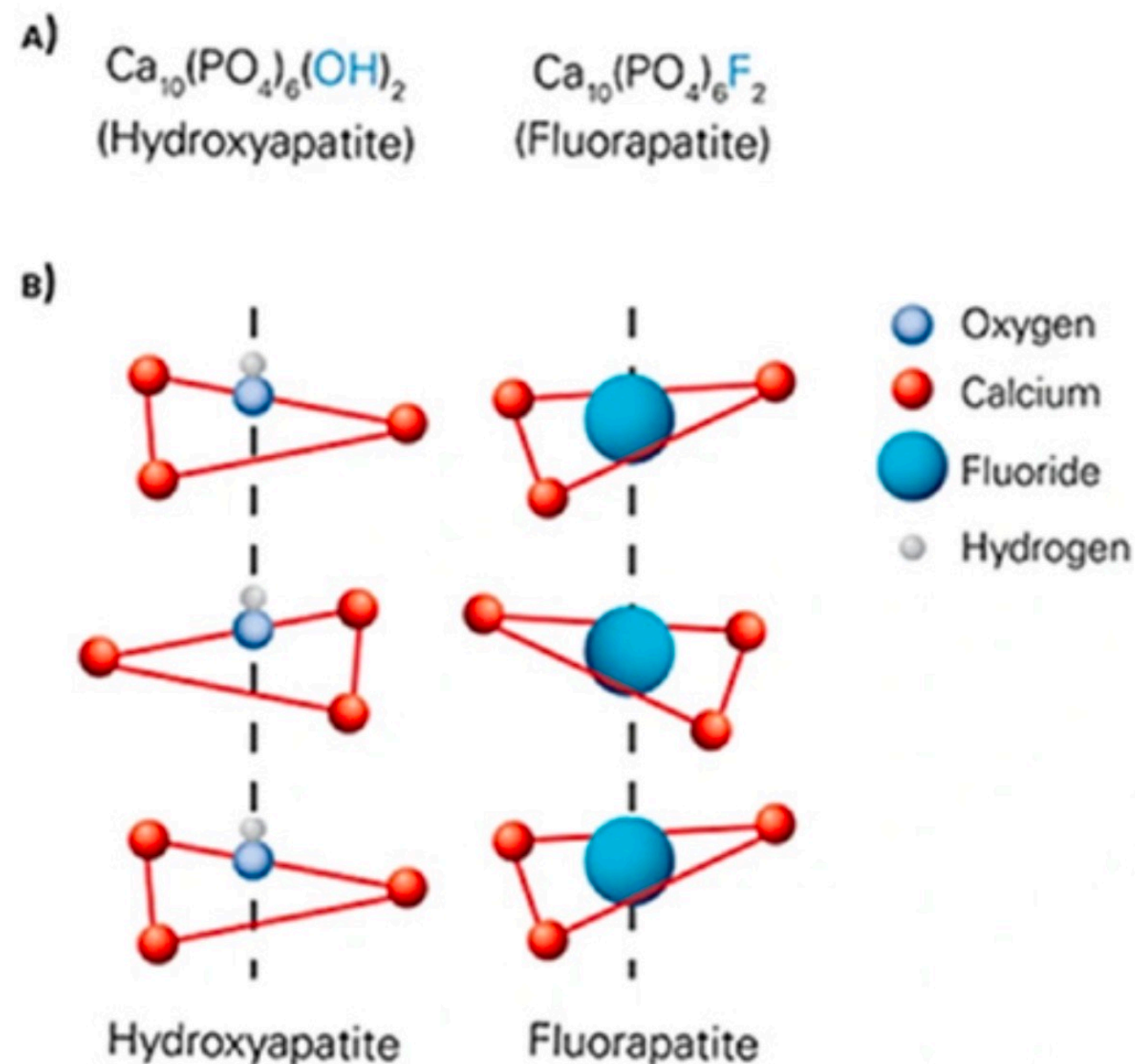
Fluoride uptake after treatment

3M™ Clinpro™ Clear Fluoride Treatment demonstrated a greater fluoride uptake on enamel.

Clinpro Clear Fluoride Treatment achieved these results in only 15 minutes.

Source: Internal. Data on file.
Laboratory testing using a model similar to FDA Method 40 for enamel fluoride uptake.

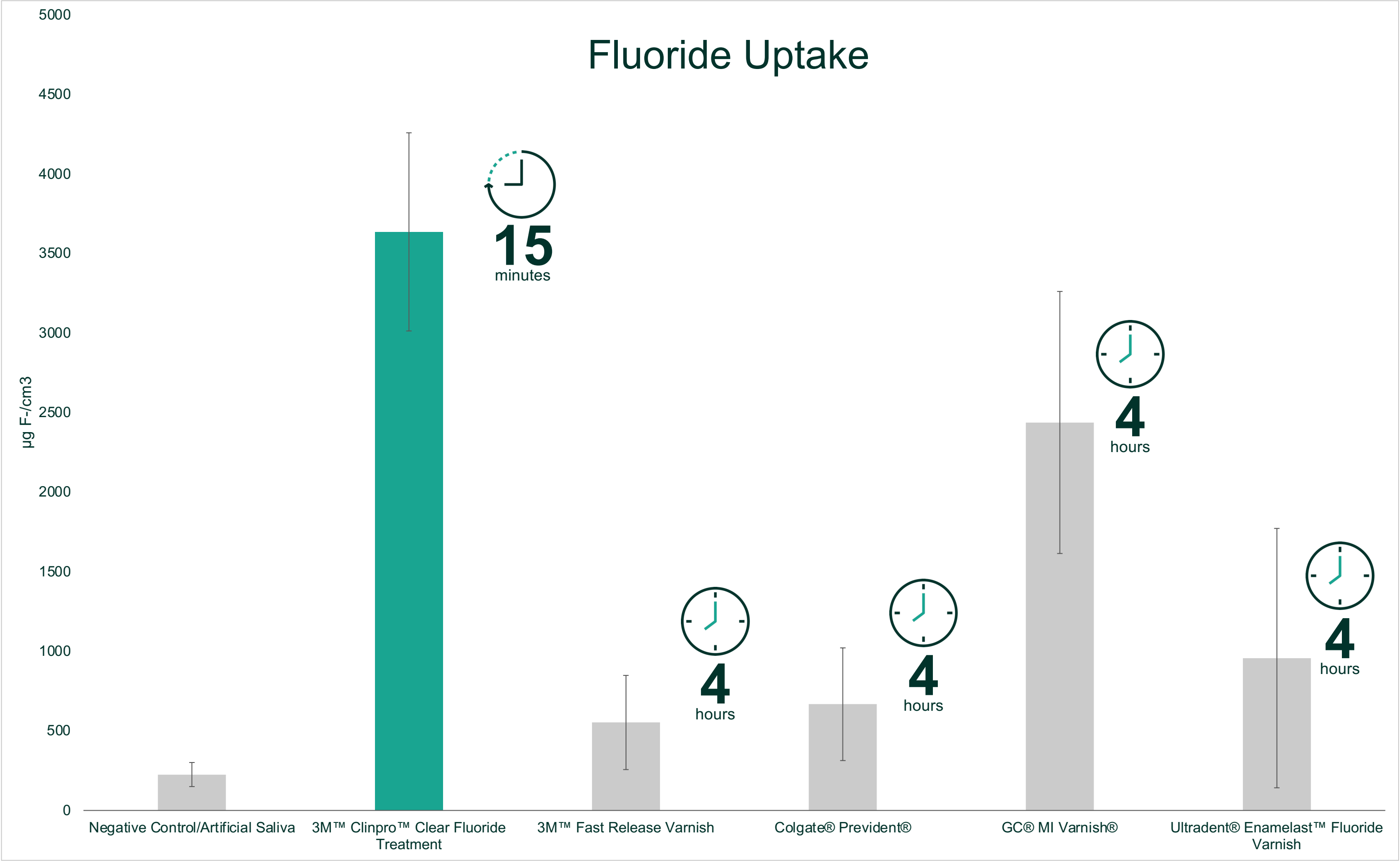
Uptake measurements were taken after elapsed treatment times specified in the instructions for use for each product.

Fluoride uptake after treatment

3M™ Clinpro™ Clear Fluoride Treatment demonstrated a greater or equal fluoride uptake on enamel.

Clinpro Clear Fluoride Treatment achieved these results in only 15 minutes.

Fluoride Uptake and Remineralization From a New Fluoride Treatment. Kuznia, M., *J Dent Res* Vol Spec Iss#104:1373

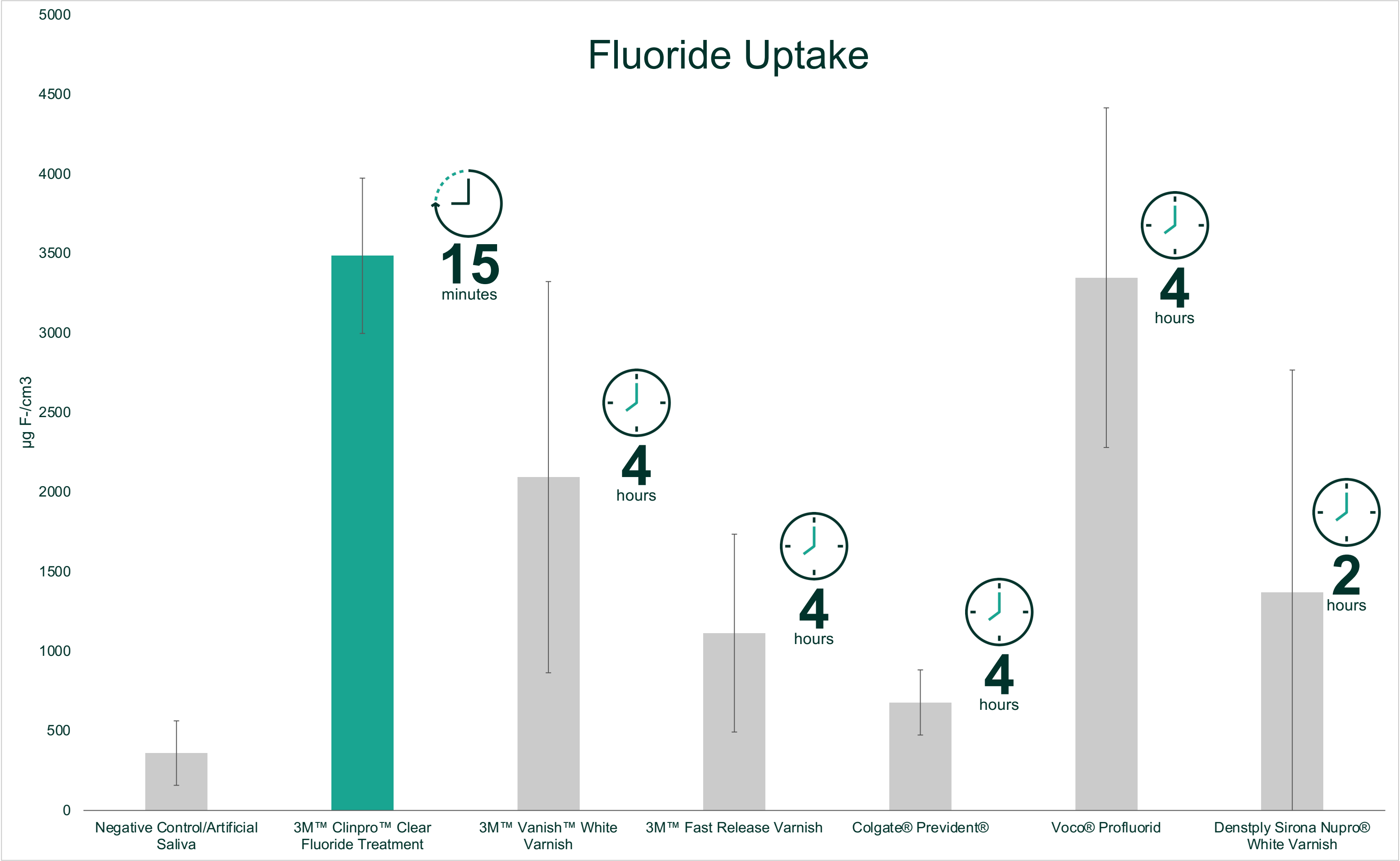
Uptake measurements were taken after elapsed treatment times specified in the instructions for use for each product.


Fluoride varnish

3M™ Clinpro™ Clear Fluoride

4
hours

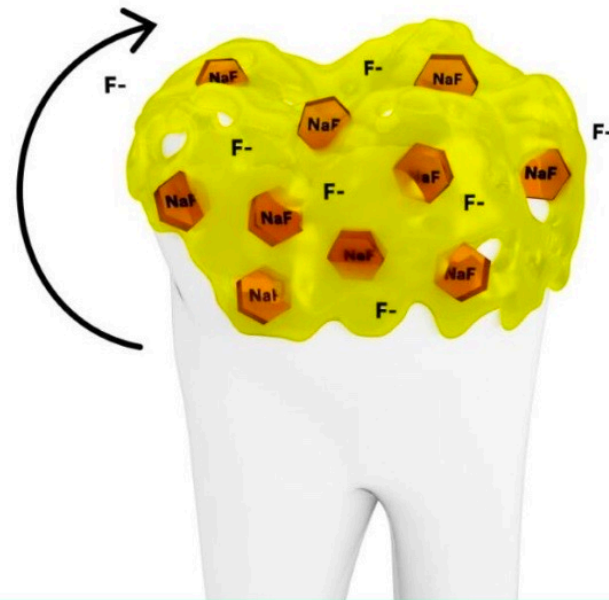Ready
release
formula

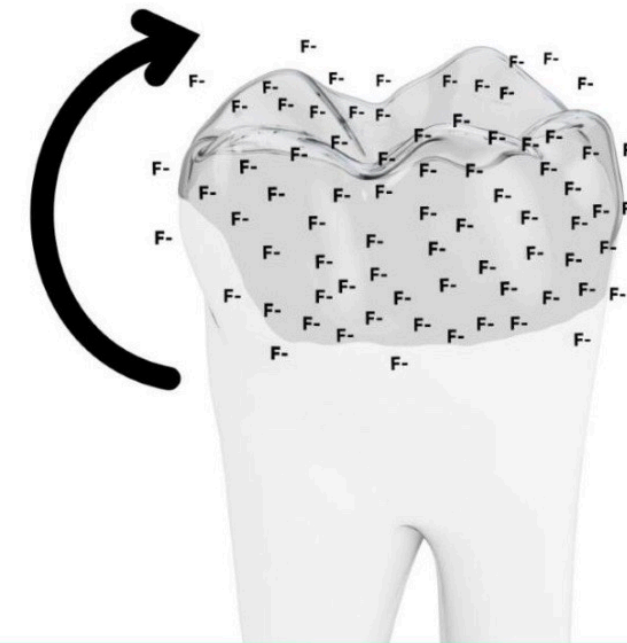15
minutes

Traditional varnishes

- Fluoride in higher concentration (22,600ppm)

But

- Must stay in contact with the tooth surface for **hours*** to allow the sodium fluoride particle to ionize. Only then can the tooth uptake the fluoride ions
- **Sodium fluoride must be dissolved from the rosin prior to uptake by the tooth**

3M™ Clinpro™ Clear Fluoride Treatment


- Fluoride in lower concentration (9,500ppm)

And

- The soluble fluoride ions are immediately available to uptake onto the tooth surface, with a treatment taking only **15 minutes**
- **Fluoride ions are soluble and ready to release**

*Fernandez et al. showed that only 258 µg F/g are free and soluble fluoride with a rosin-based varnish, the rest is insoluble. That means a longer retention time (4 hours) is needed to allow NaF particles initially insoluble in the varnish matrix to release for uptake on the tooth surface. Brazilian Dental Journal (2014) 25(2): 160-164

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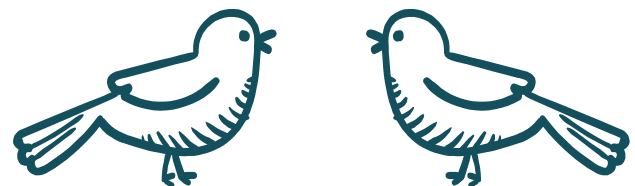Fluoride varnish

3M™ Clinpro™ Clear Fluoride

Water-based technology with ionic minerals

- 3M™ Clinpro™ Clear Fluoride Treatment attaches to the tooth using an ion exchange reaction
- Minerals are fully ionic in solution and immediately available, removing the need to wait for sticky resins to release them at the tooth's surface
- Contains fluoride, calcium, and phosphate

Benefit: Minimal contact time
for effective performance.


Fluoride varnish

Fluoride concentrations in professionally applied treatment agents

Topical fluoride agents	PPM fluoride	Percent fluoride	Percent sodium fluoride
2.26% fluoride varnish ²	22,600	2.26	5.0
Acidulated phosphate fluoride (APF) with phosphoric acid ³	12,300	1.23	2.7
3M™ Clinpro™ Clear Fluoride Treatment	9,500	0.95	2.1 ¹
Prescription gel or paste	5,000	0.50	1.1
Prescription fluoride mouth rinse	900	0.09	0.2

¹ Double concentration of prescription gels.
² Available fluoride in varnish is much lower because fluoride is solid rather than solution.
³ Acidic products will also dissolve mineral and can etch esthetic dental work and cause pitting. Can also burn and irritate.

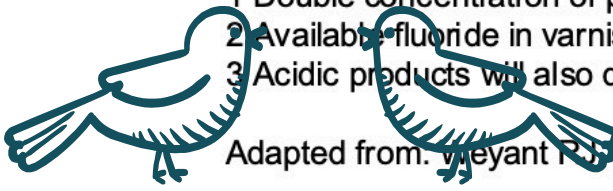Fluoride varnish



3M™ Clinpro™ Clear Fluoride

Fluoride uptake equivalent to traditional varnishes in less time

Samples were tested following manufacturer's instructions for use on treatment time.


The amount of fluoride deposited on teeth using 3M™ Clinpro™ Clear Fluoride Treatment shows no statistical difference to traditional fluoride varnishes.

Clinpro Clear Fluoride Treatment achieves these results in only 15 minutes.

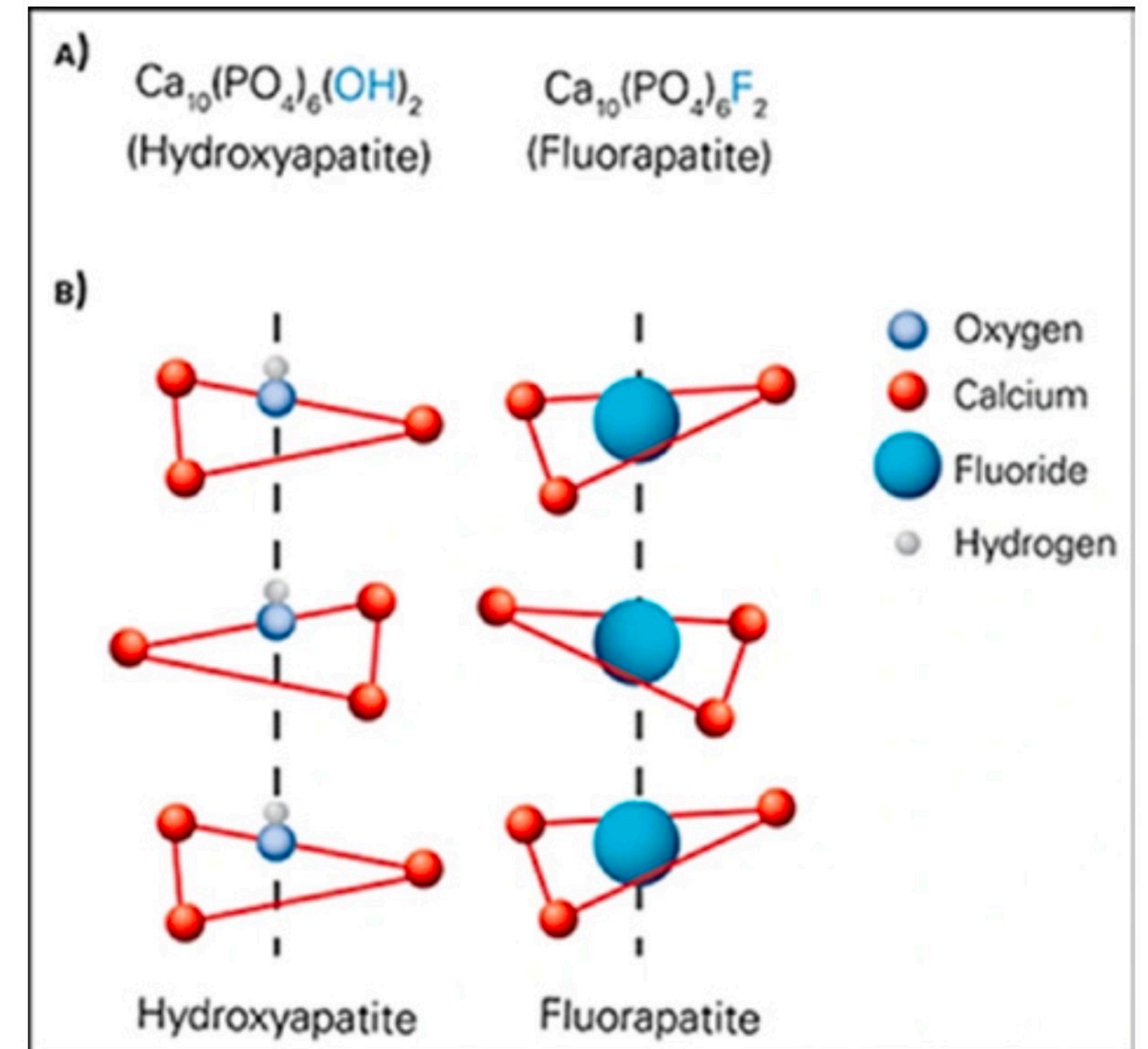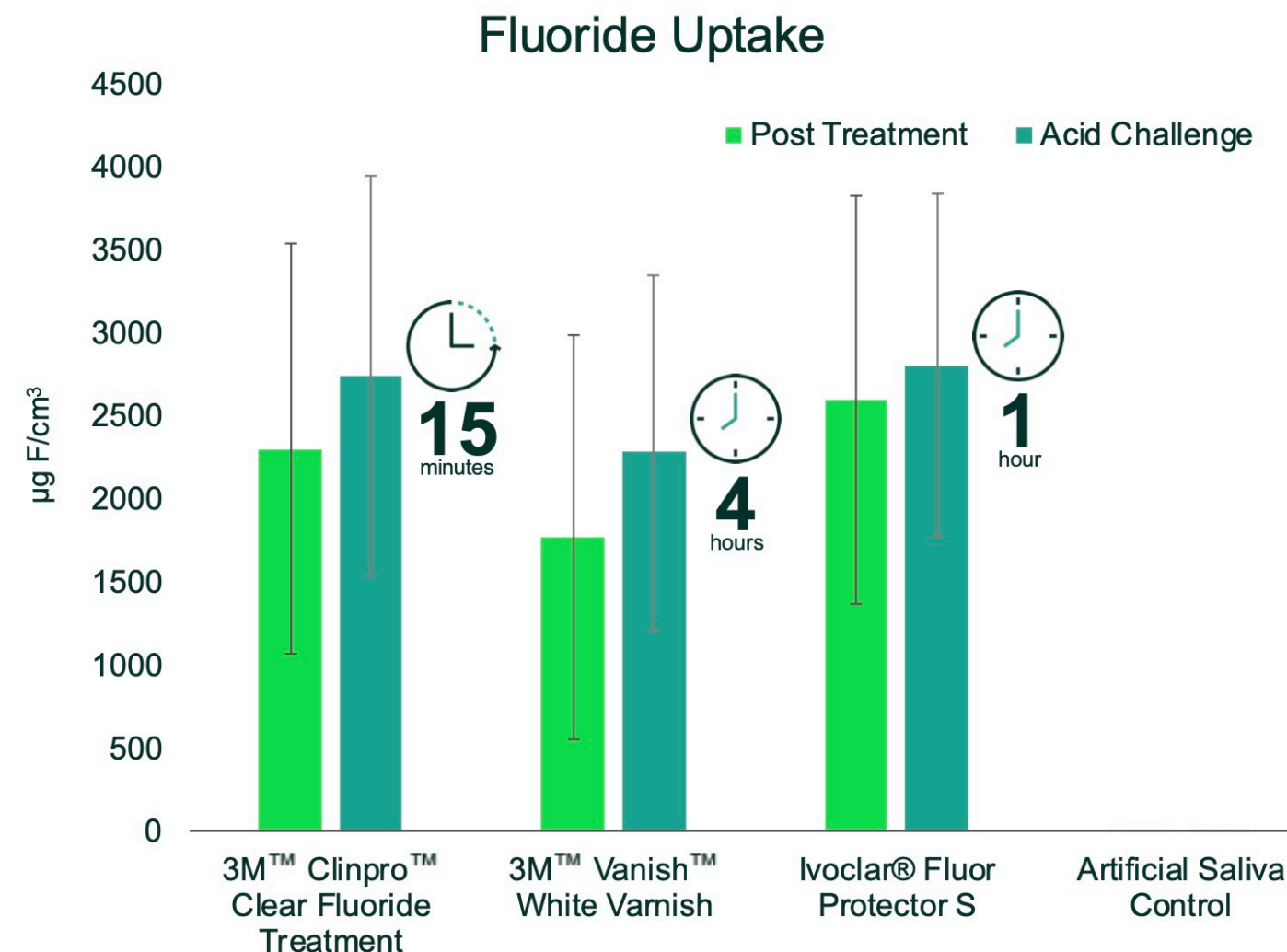Source: Internal. Data on file.
Laboratory testing using a model similar to FDA Method 40 for enamel fluoride uptake

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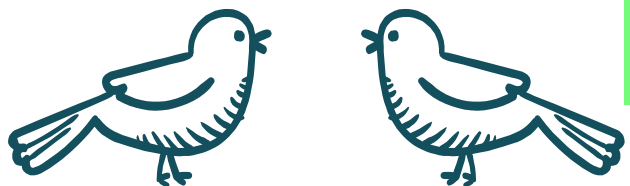The relevance of soluble fluoride


3M™ Clinpro™ Clear Fluoride

In vitro study at the University of Michigan comparing the concentration of total and soluble fluoride among different fluoride treatments.

For total fluoride:

- 3M™ Clinpro™ Clear Fluoride Treatment showed a total fluoride concentration around 9,500 µg F/g
- All other varnishes had a total fluoride concentration close to 22,600 µg F/g

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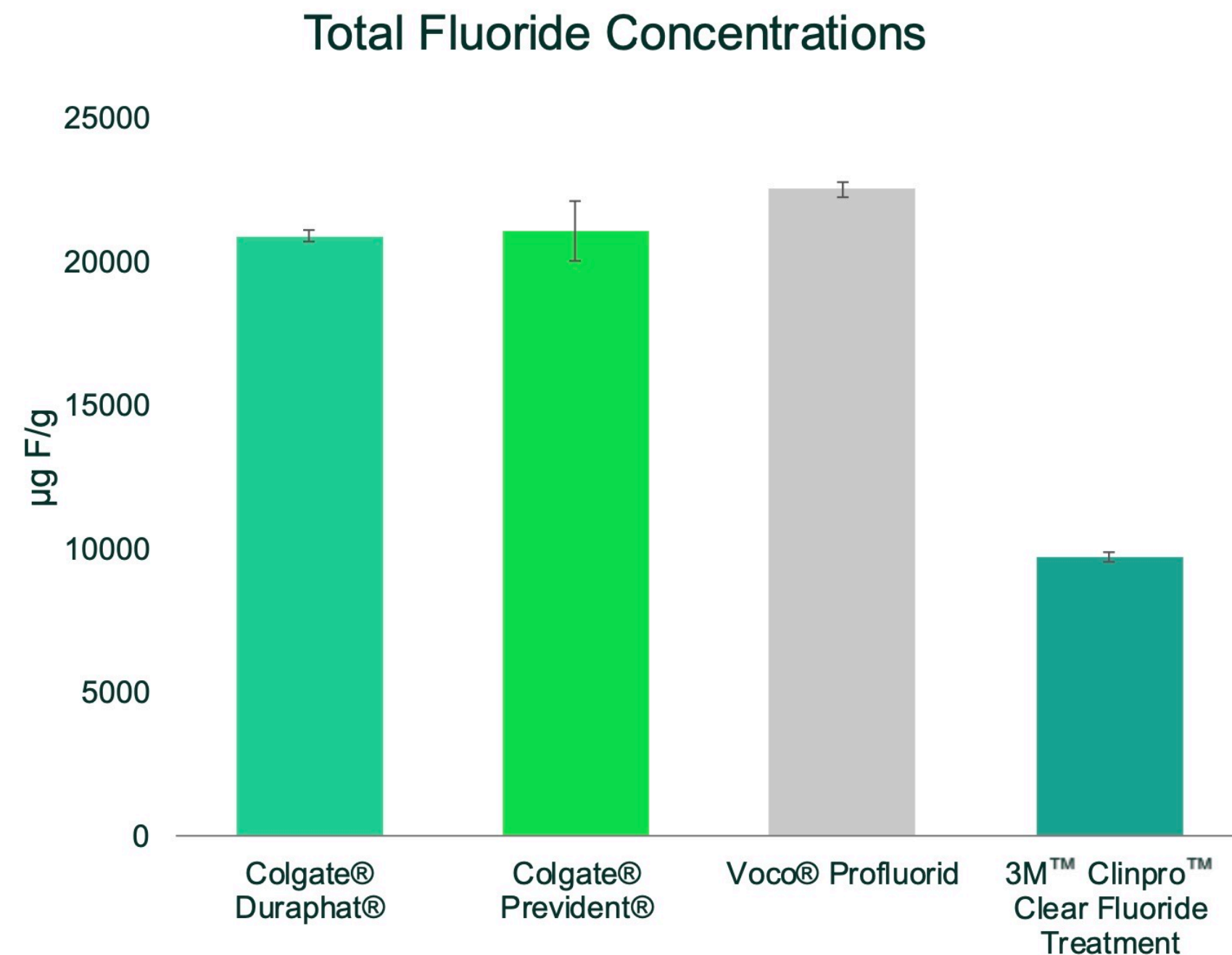
The relevance of soluble fluoride

3M™ Clinpro™ Clear Fluoride

The study then compares the amount of soluble fluoride in each product

- 3M™ Clinpro™ Clear Fluoride Treatment still measured a concentration of fluoride around 9,500 µg F/g, showing 100% of its fluoride is in a soluble form
- In the fluoride varnish groups, measured fluoride concentrations for soluble fluoride were a low percentage of the total measured fluoride

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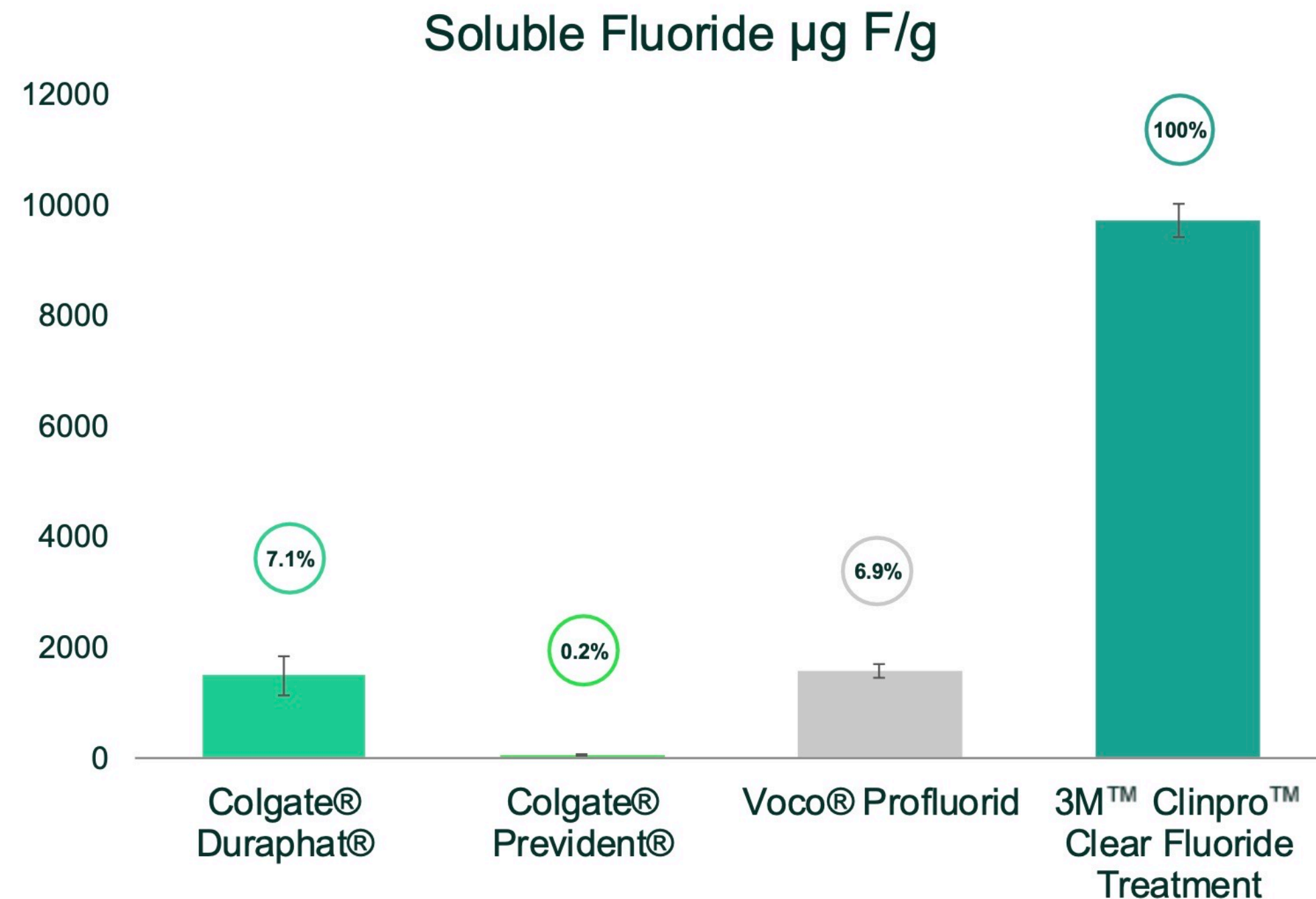Source: University of


Fluoride varnish 3M™ Clinpro™ Clear Fluoride

The relevance of soluble fluoride

Even with a lower overall fluoride concentration, 3M™ Clinpro™ Clear Fluoride Treatment provides a greater quantity of available fluoride to the tooth compared to rosin-based varnishes.

9,500 µg F/g goes in, and 9,500 µg F/g is released in solution

Clinpro Clear Fluoride Treatment is a ready release formula

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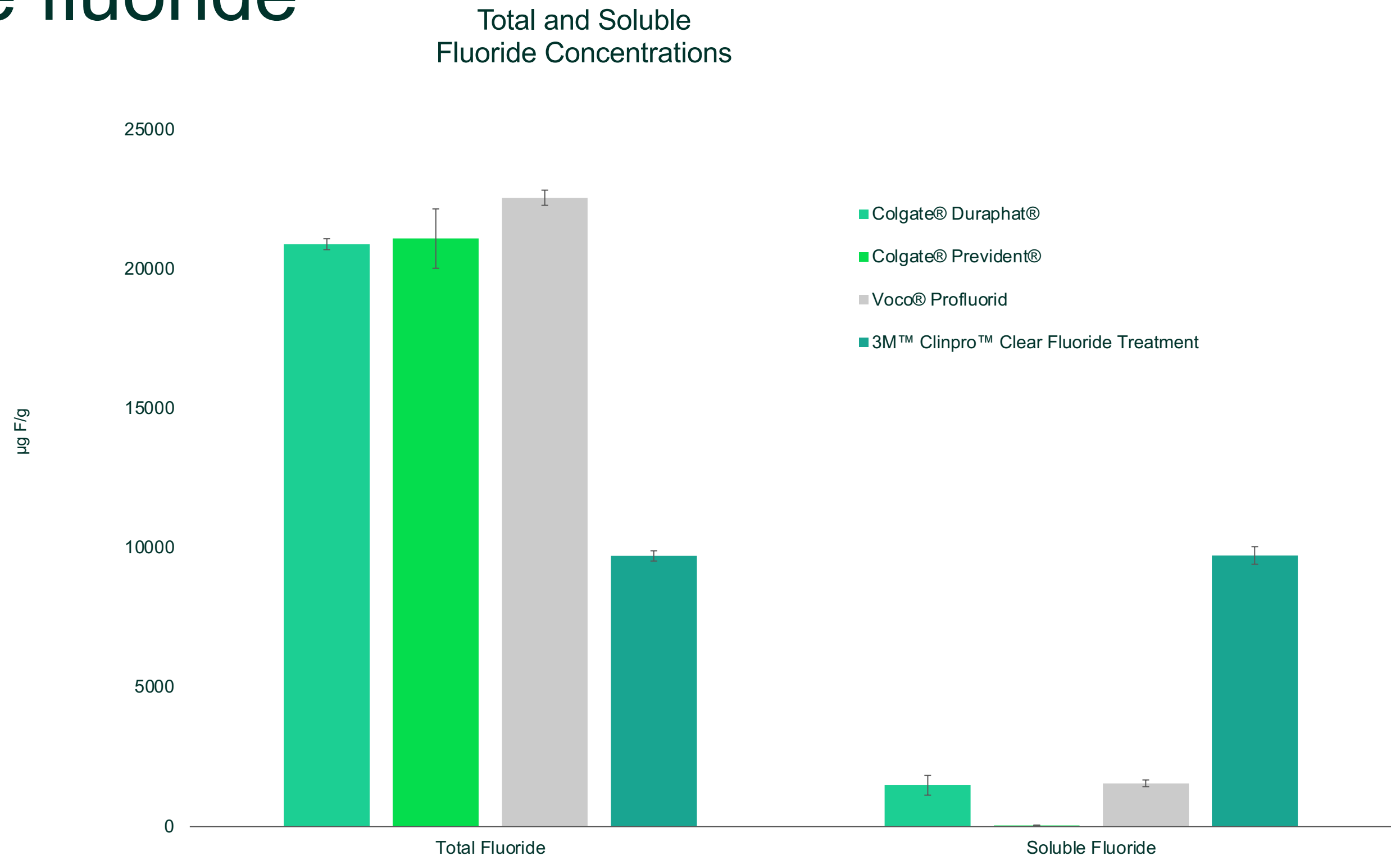
The relevance of soluble fluoride

Fluoride deposits

- Loosely bound fluoride reservoirs were measured
- Once dissolved, these reservoirs can increase the fluoride concentration in the microenvironment on the surface of teeth
- In demineralized enamel, 3M™ Clinpro™ Clear Fluoride Treatment resulted in the highest deposited fluoride compared with all other products

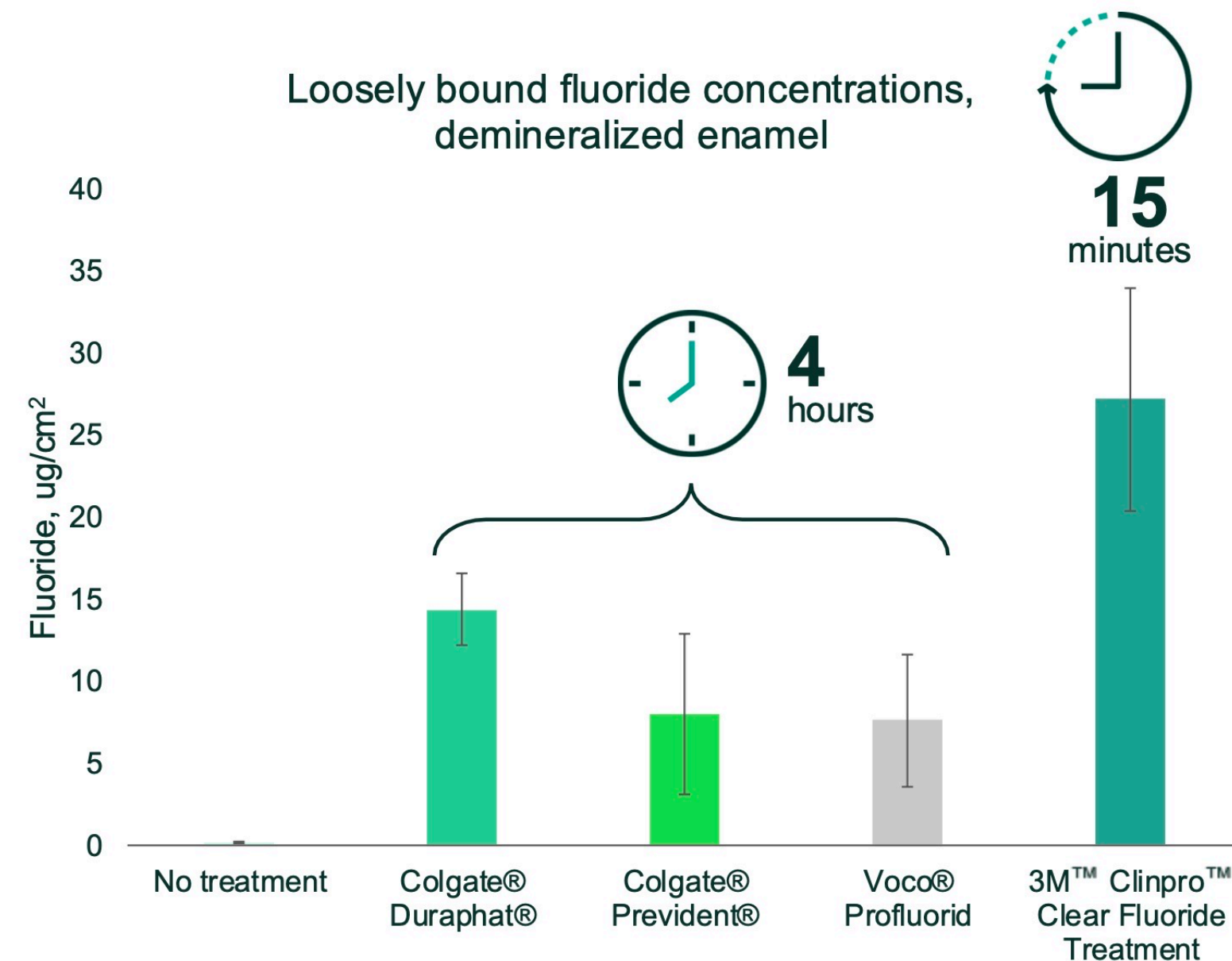Source: University of Michigan Report, Dr. Livia Tenuta. Data on file.

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Fluoride varnish

3M™ Clinpro™ Clear Fluoride

Table 1. FLUORIDE IN THE TESTED PRODUCTS, RELEASE TO ARTIFICIAL SALIVA, AND REACTIVITY TO SOUND AND DEMINERALIZED ENAMEL (AVERAGE±SD)*

Product	Total fluoride (ppm F, N=3)	Fluoride available in the soluble fraction (ppm F, N=3)	Fluoride release to artificial saliva (ppm F, N=5)	Fluoride reactivity**	
				Sound enamel (mg F/cm ² , N=12)	Demineralized enamel (mg F/cm ² , N=8)
Negative control†	-	-	0.031±0.002‡	0.17±0.08	0.22±0.05
Clinpro™ Clear	9,713±184	9,728±316	5.844±1.209	5.06±2.32	27.2±6.77
Duraphat®	20,892±195	1,485±353	0.043±0.005‡	3.64±1.63	14.4±2.21
Prevident®	21,093±1,065	50±9	0.349±0.149	2.76±0.76	8.00±4.92
Profluorid®	22,558±272	1,558±120	0.078±0.011	0.98±0.54	7.63±4.03

* Abbreviations in table: SD=standard deviation; ppm F=parts per million Fluoride.

** Alkali-soluble fluoride concentration.

† Negative control=enamel slabs which were not treated, but were cycled through artificial saliva as the other groups and had their enamel fluoride concentration determined after that.

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3M™ Clinpro™ Clear Fluoride

Less is more

Teeth receive the 2.1% dose of sodium fluoride

Traditional 5% NaF varnishes

- The fluoride ions are trapped in suspension until the rosin dissolves
- Like blending salt into honey to create a sticky mixture. The mixture will taste sweet until the honey dissolves, releasing the salt to your tastebuds
- Time must elapse to release the salt

3M™ Clinpro™ Clear Fluoride Treatment

- The fluoride ions are freely available in a water-based solution
- Like dissolving salt in a glass of water. All the salt dissolves and the solution tastes salty right away
- The salt is fully soluble and ready to release

Fluoride varnish

3M™ Clinpro™ Clear Fluoride

A whole new way to deliver fluoride treatments

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 Smooth, water-based

 2.1% sodium fluoride

 Stabilized calcium phosphate technology

 15-minute contact time

 Unique L-Pop delivery system

Fluoride varnish


3M™ Clinpro™ Clear Fluoride

Benefits for patients and professionals

การปรับปรุงของ 3M™ Clinpro™ Clear Fluoride จาก Fluoride เดิม			
① Improved patient experience	② Formulation	③ Efficient fluoride uptake	④ Easy application
• Short contact time • Pleasant taste, smell, and mouthfeel • Non-sticky • More likely to say “yes” to treatment	• Water-based • Stabilized calcium phosphate technology • Neutral pH • Rosin-free • Free of common allergens*	• 2.1% vs. 5% NaF • Ready release, high uptake • Equivalent performance to traditional treatments • 15-minute contact time	• Unique L-Pop delivery system • No mixing • Just press and go

*See allergen disclosure at 3M.com/ClinproClear

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



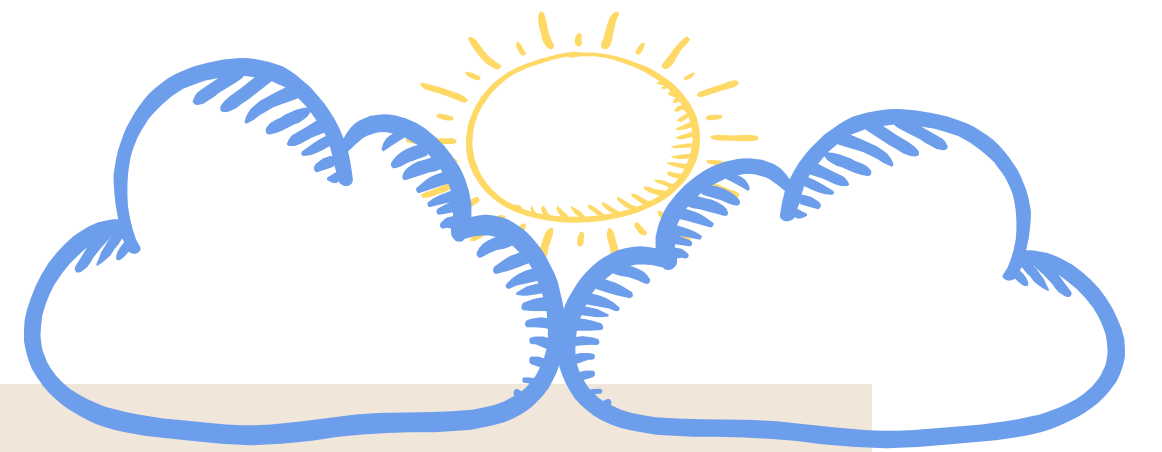
- Active white lesion
- Prone to caries/difficult to access
- Occlusal surface of newly erupted permanent molar

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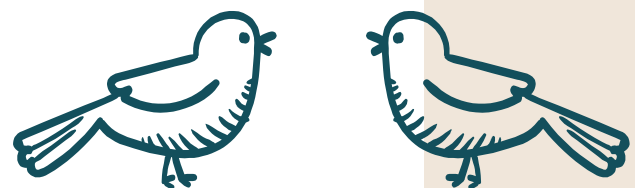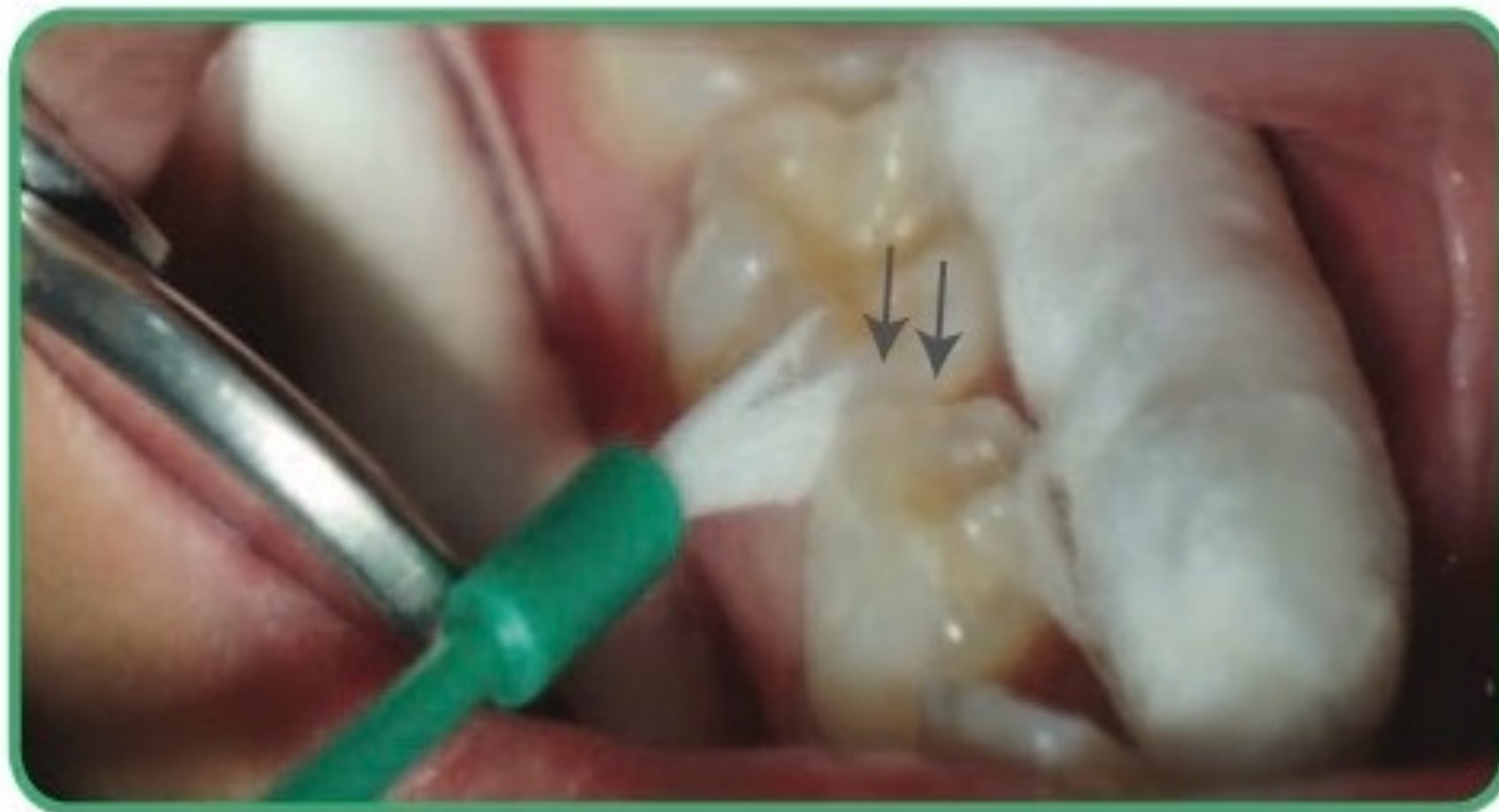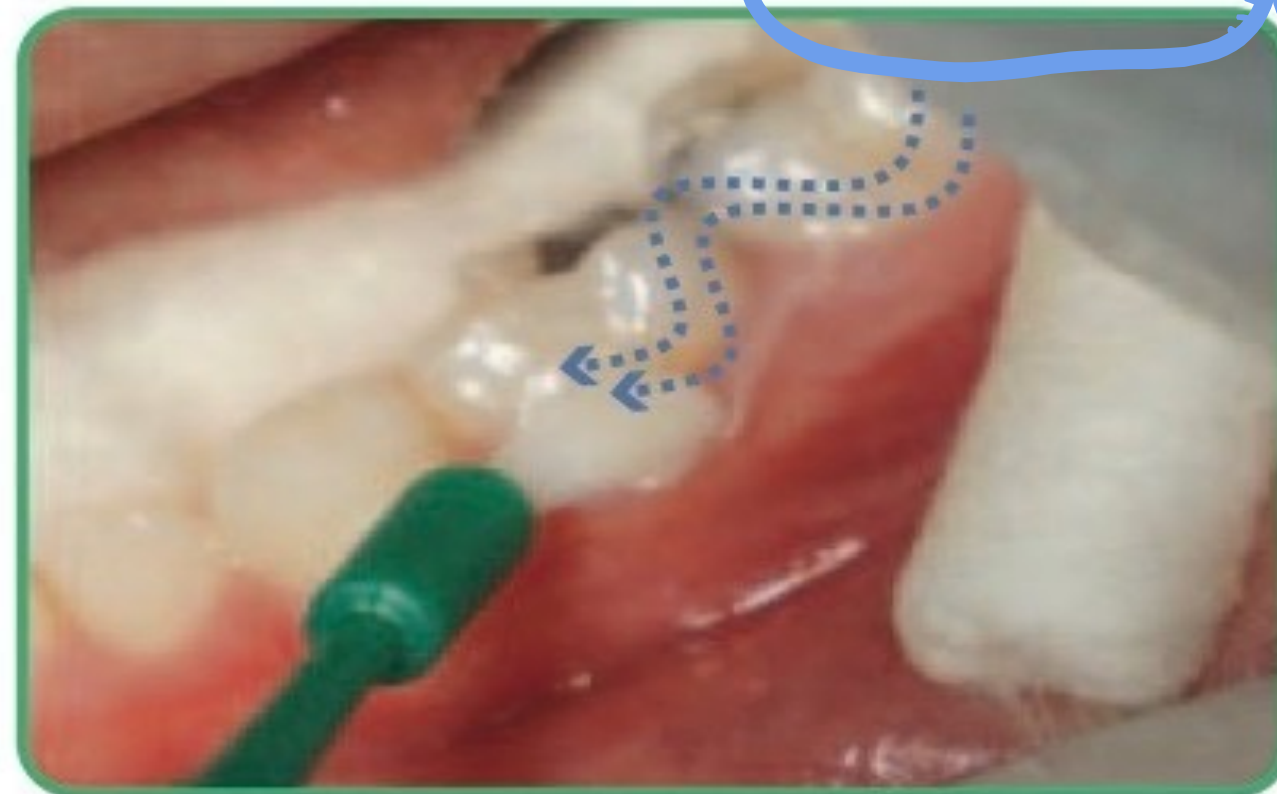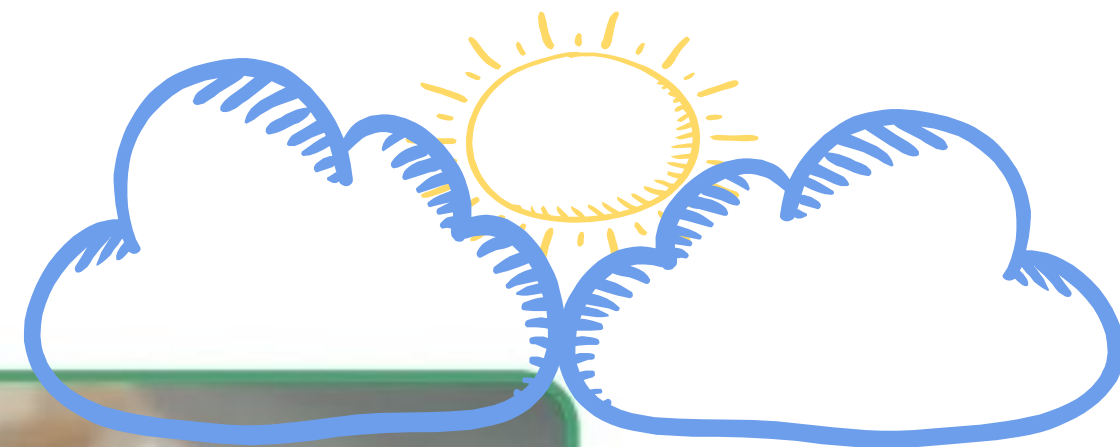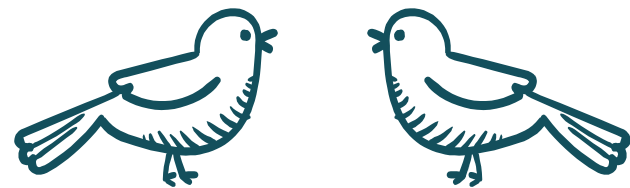
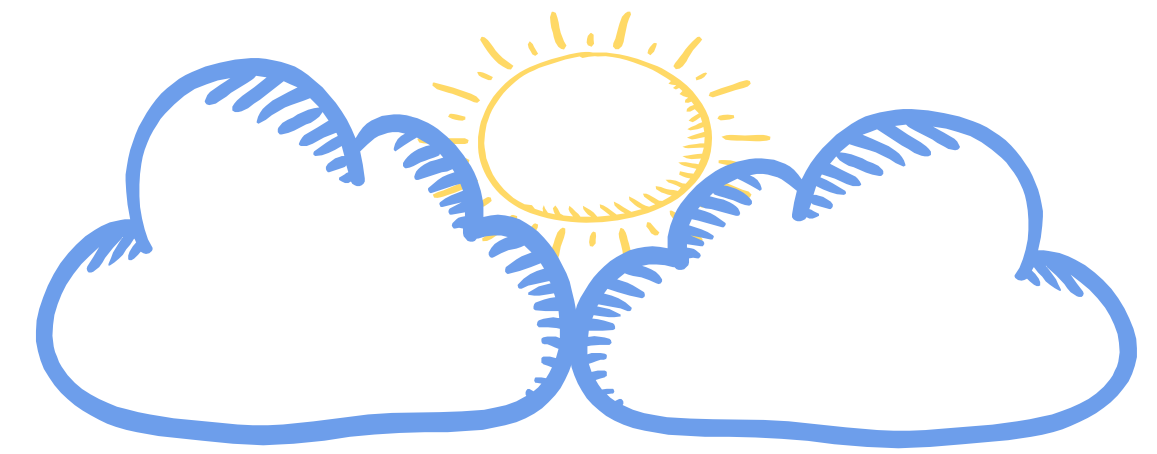
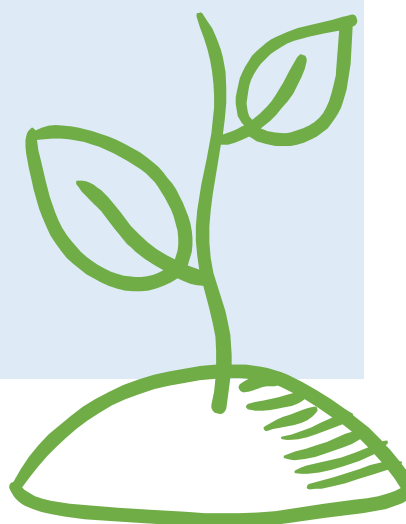
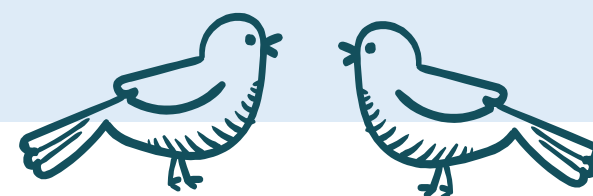
Post application instructions

Note : No fluoride supplement
and topical high fluoride products

	DAT (Thailand) 2565
Avoid drinking	-
Avoid eating soft food	-
Avoid eating hard food	2 - 4 hr
Avoid brushing/flossing	All day

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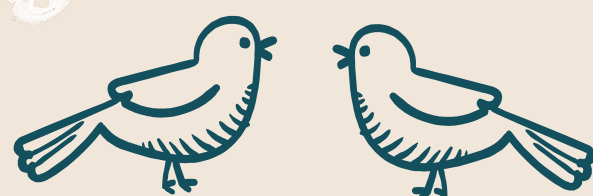
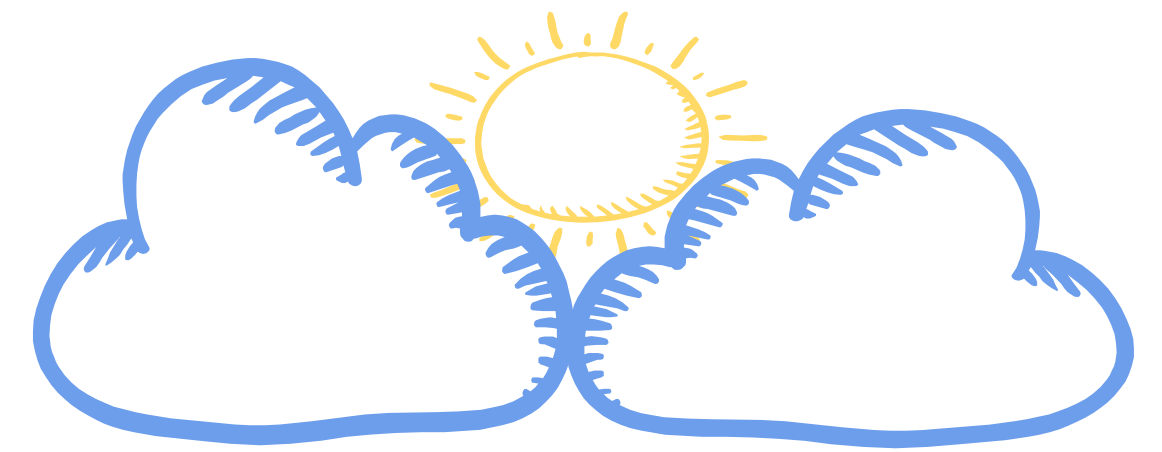
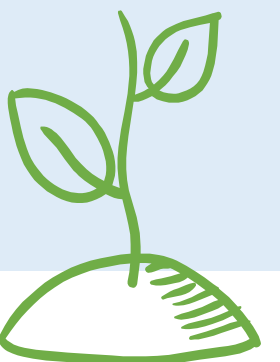
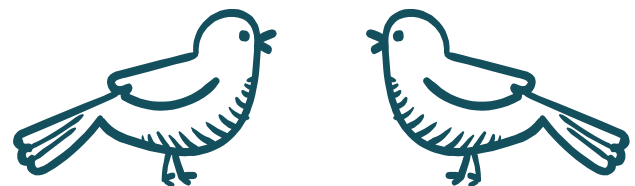
3M™ Clinpro™ Clear Fluoride

Relief for sensitive teeth

When dentin tubules are exposed, stimuli can cause fluid flow and pressure on nerves, causing pain or sensitivity. Dentin tubule exposure can be caused by:

- Gingival recession
- Erosive wear
- Periodontal treatment

Open dentin tubules

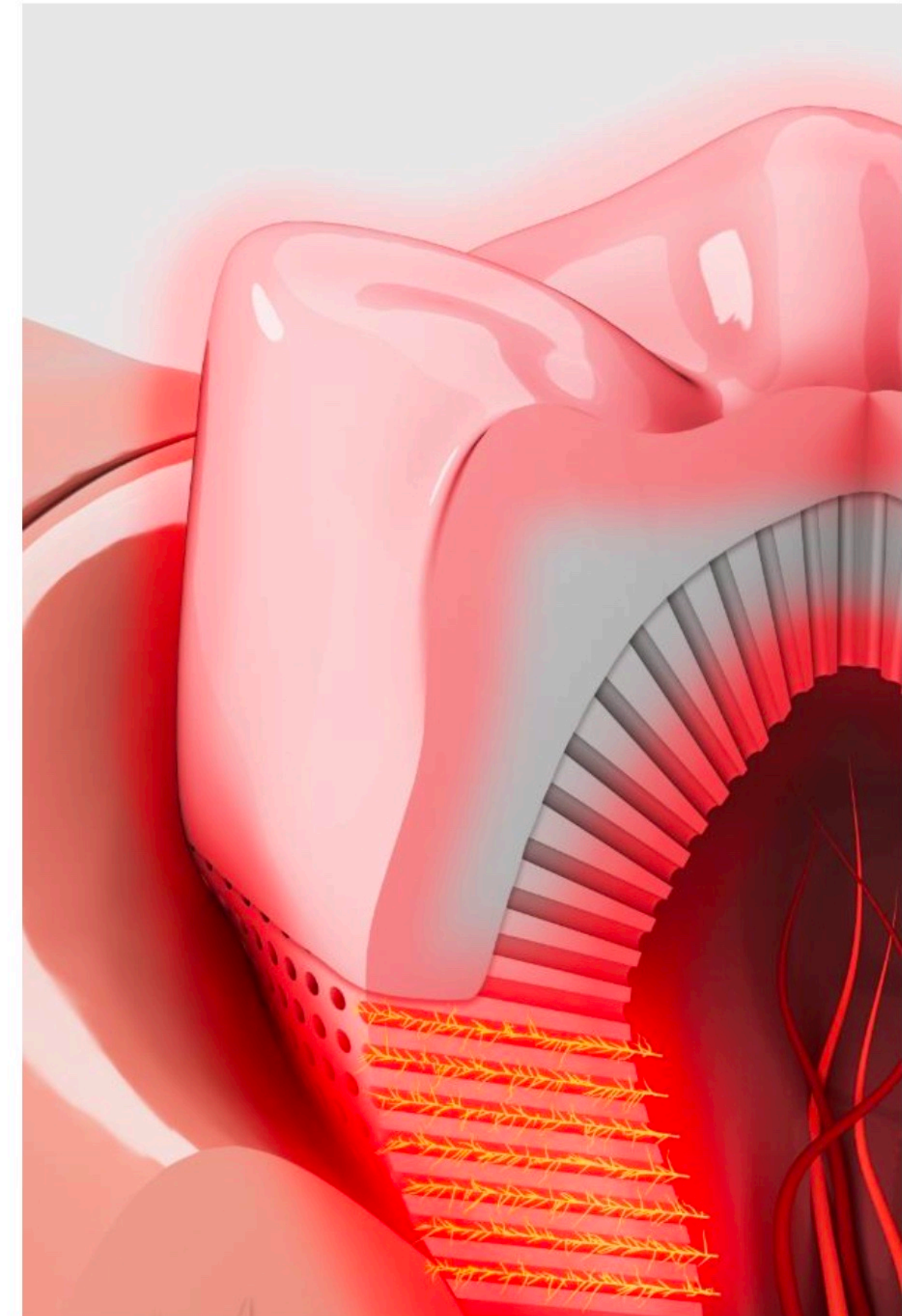Fluoride varnish

3M™ Clinpro™ Clear Fluoride



Relief for sensitive teeth

Applying 3M™ Clinpro™ Clear Fluoride Treatment creates a physical barrier that reduces hypersensitivity.

After application, Clinpro Clear Fluoride Treatment forms mineral deposits which block dentin tubules and reduce dentin fluid flow.

This reduction in dentin fluid flow reduces hypersensitivity.

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
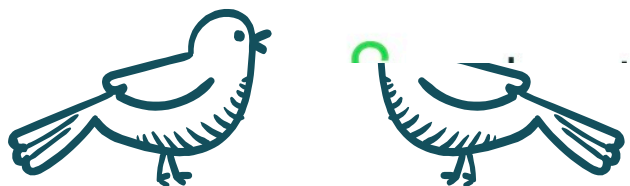
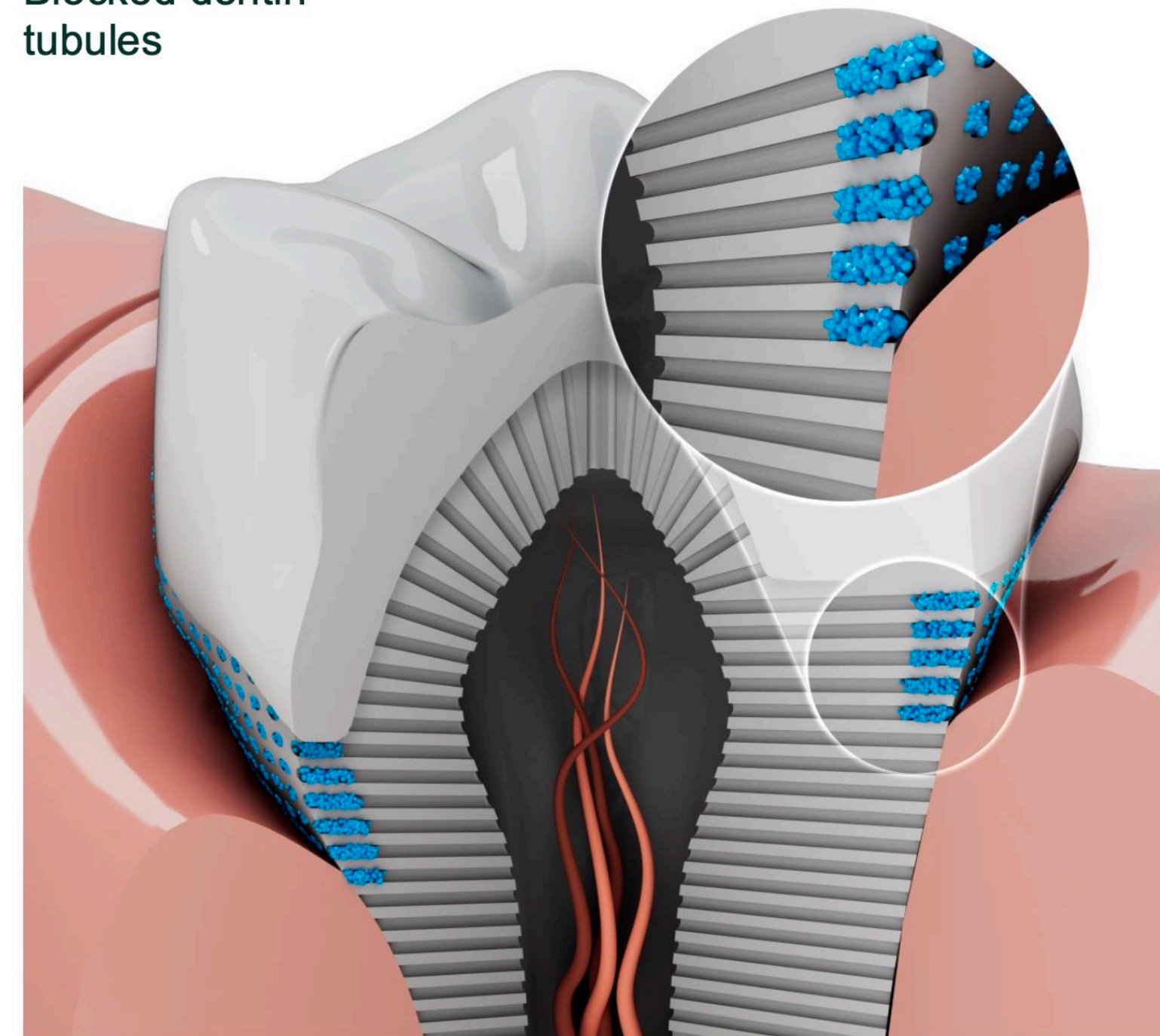
“TIPS AND TRICKS WITH
CLINPRO CLEAR & SEALANT”


- 6 AUG 2025 -

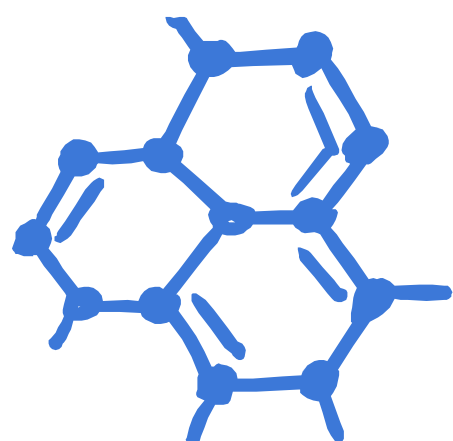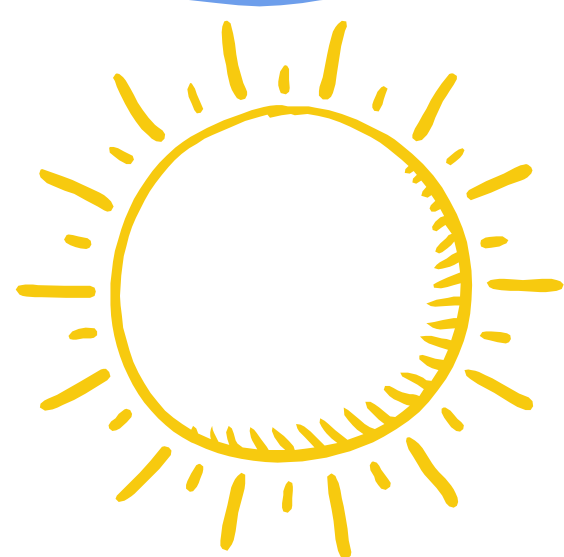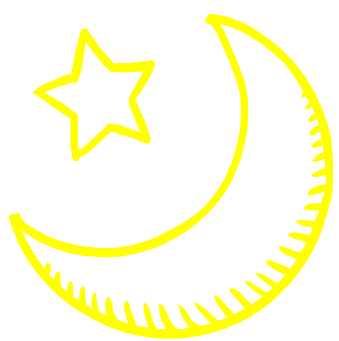Thanks you

δ solventum

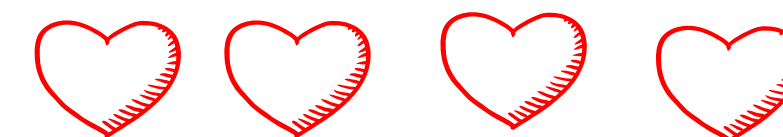Reference

Scientific Meeting Publications				
Publication Year	Publisher	Format	Title	Link
2019	IADR	Abstract	Fluoride-Uptake, remineralization, and Micro-CT characterization of demineralization resistance of a novel in-office fluoride treatment	2019 IADR Abstract: Fluoride Uptake, Remineralization, and Micro-CT
2019	IADR	Abstract	Adhesion Bond-Strength Measurements After a Novel in-Office Fluoride Treatment	2019 IADR: Adhesion Bond-Strength Measurements After a Novel in-Office Fluoride Treatment
2022	IADR	Abstract	Cumulative Ion Release Comparison of Two Popular Fluoride Varnishes	2022 IADR Abstract: Cumulative Ion Release Comparison
2023	AADOCR	Abstract	ISO Fluoride Release Potential Study of Different Fluoride-Containing Products	2023 AADOC: ISO Fluoride release potential
2023	IADR/LAR	Abstract	Migration Potential of an Experimental Aqueous in-Office Fluoride Treatment	2023 IADR/LAR: Migration Potential of an Experimental Aqueous in-Office Fluoride Treatment
2023	AADOCR	Abstract	Remineralization and Enamel Fluoride Uptake From a New Fluoride Treatment	2023 AADOCR: Remineralization and Enamel Fluoride Uptake
2024	IADR	Abstract	Fluoride Uptake and Remineralization From a New Fluoride Treatment	2024 IADR: Fluoride Uptake and Remineralization
2024	IADR	Abstract	Retention-time and Migration-Potential of a New In-Office Fluoride Treatment	2024 IADR: Retention-time and migration
2024	IADR	Abstract	Fluoride Uptake and Remineralization From a New Fluoride Treatment	2024 IADR Abstract: Fluoride Uptake and Remineralization
2024	IADR	Abstract	Ion Release Potential of Different Fluoride-Containing Products	2024 IADR Abstract: Ion release potential
2024	IADR	Abstract	Dental Professional Preference Among Three in-Office Fluoride Treatments	2024 IADR Abstract: Dental Professional Preference
2024	IADR	Abstract	Patient Preferences Among Three in-Office Fluoride Treatments	2024 IADR Abstract: Patient Preferences
2025	AADOCR	Abstract	Fluoride Treatment Performance of Hydrophilic/Hydrophobic Compositions	
2025	AADOCR	Abstract	Effect of New Fluoride Treatment on Orthodontic Bracket Bond Strength	
2025	AADOCR	Abstract	Fluoride Uptake and Remineralization from a New Fluoride Treatment	


External/3rd Party				
Publication Year	Publisher	Format	Title	Link
2025	Pediatric Dentistry	Journal	Enamel fluoride reactivity of professional fluoride products is not correlated with their total fluoride content	https://pubmed.ncbi.nlm.nih.gov/40533917/
2025	Dental Materials ADM Supplements 2023	Abstract	157 - Characterization of Remineralized Artificial Enamel Lesion using SEM	https://doi.org/10.1016/j.dental.2025.01.170

Fluoride varnish


3M™ Clinpro™ Clear Fluoride

Relief for sensitive teeth


Uncoated bovine dentin

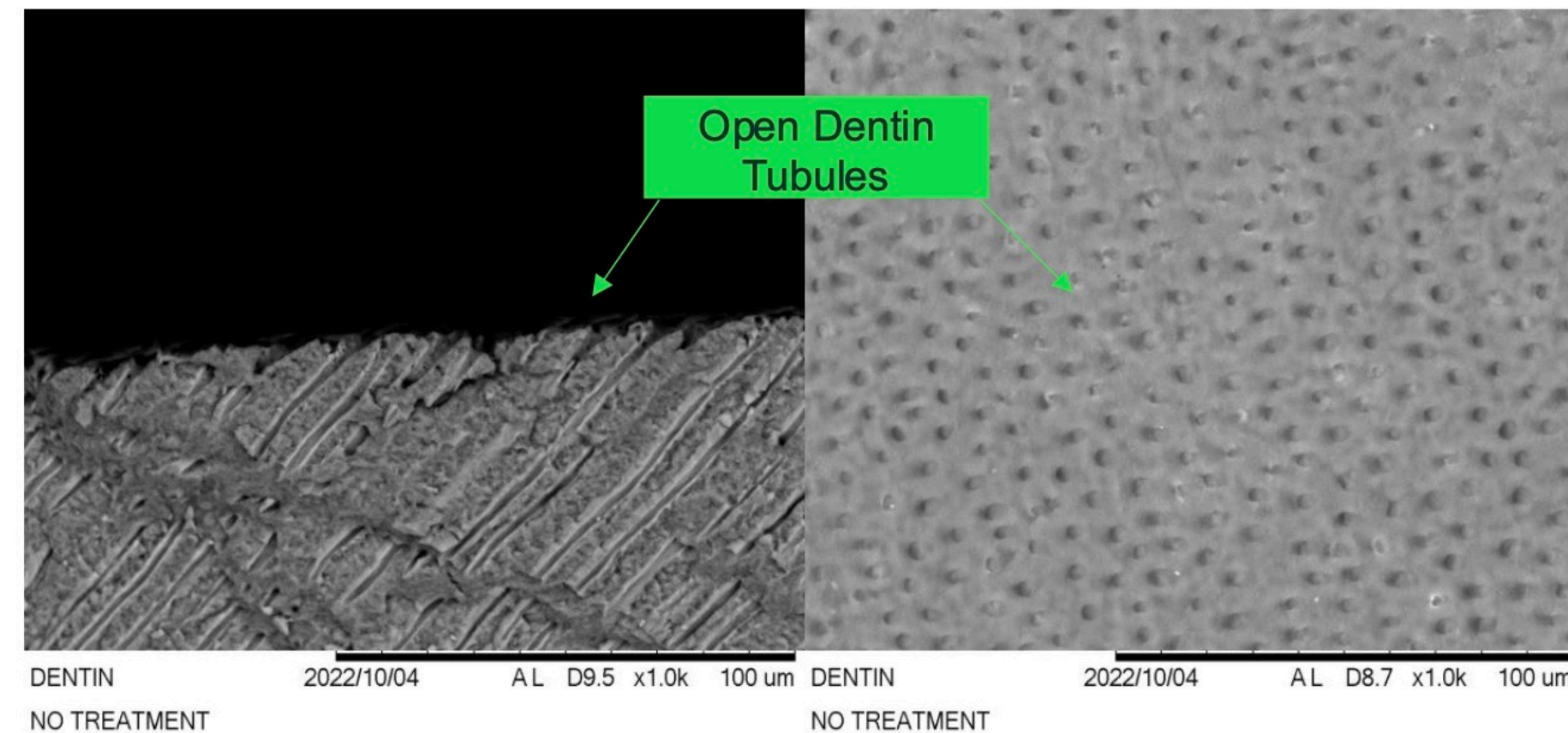
Bovine dentin coated with
3M™ Clinpro™ Clear Fluoride Treatment
(In 15 minutes artificial saliva)

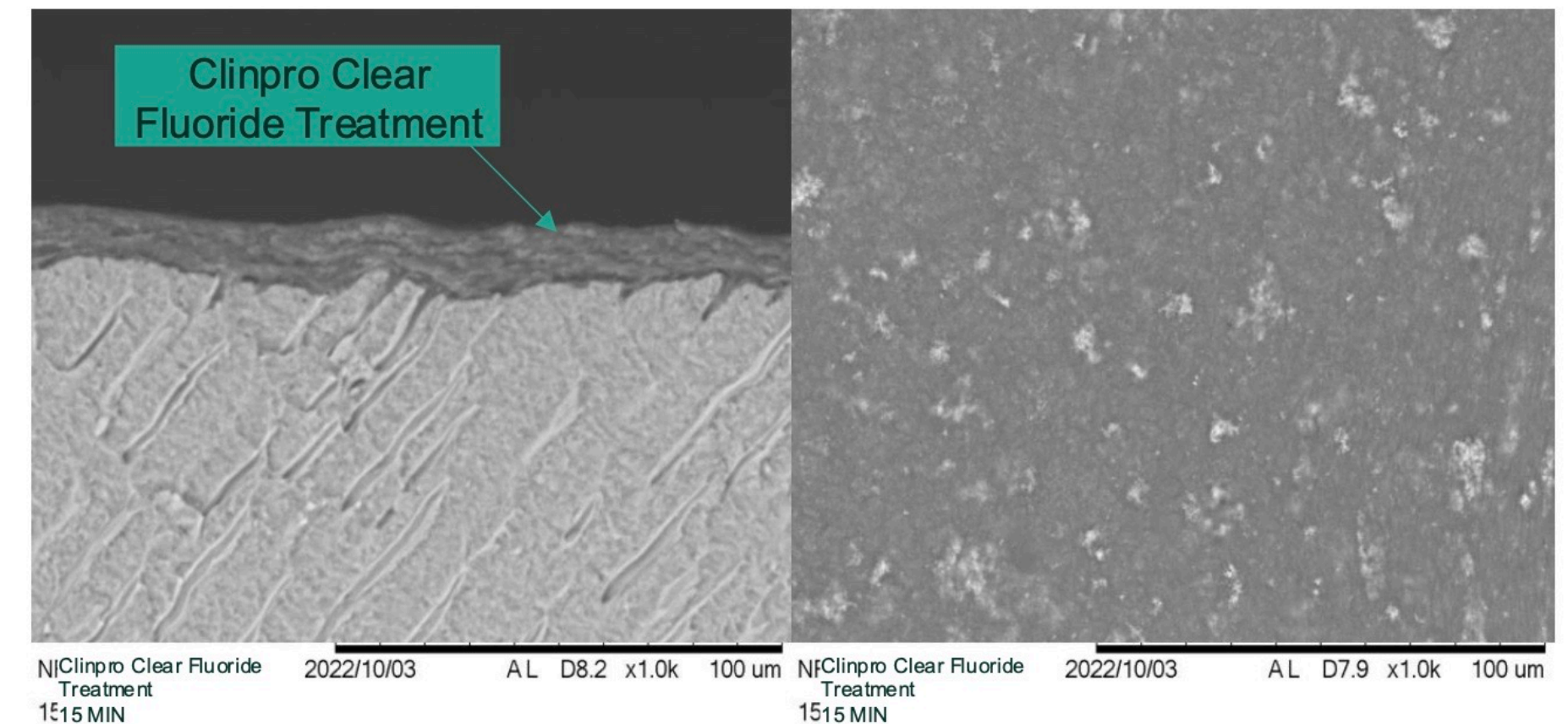Ready release formula

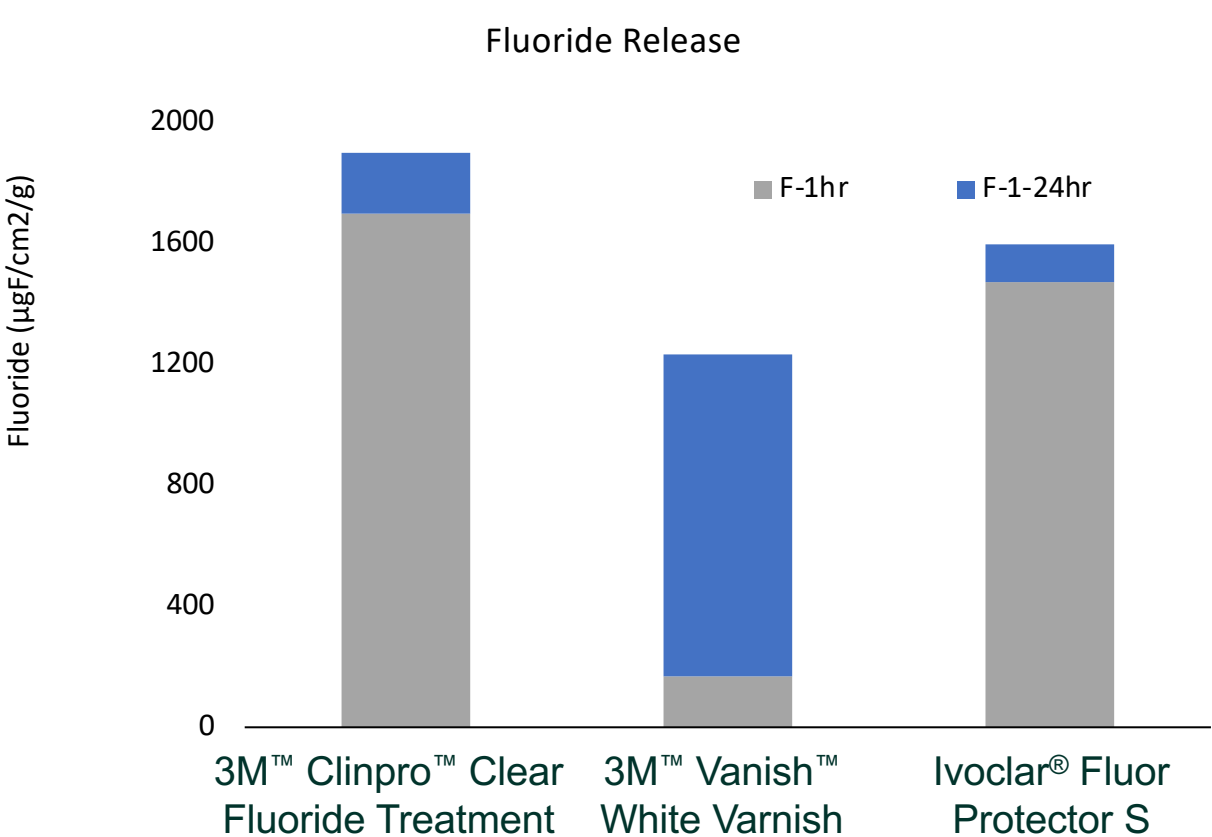
- The cumulative fluoride release of Clinpro Clear Fluoride Treatment at 24 hours is statistically higher than Vanish White Varnish, and statistically not different from Ivoclar Fluor Protector S

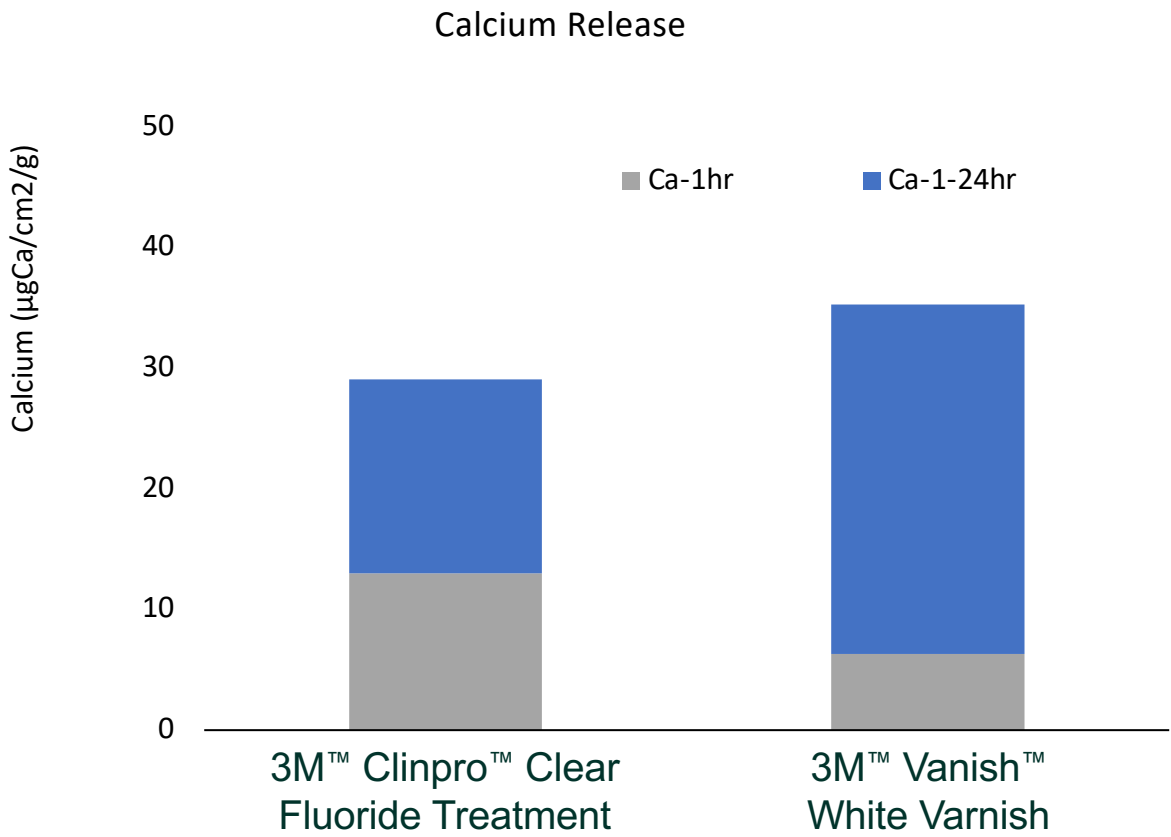- The cumulative calcium release of Clinpro Clear Fluoride Treatment is statistically not different than Vanish White Varnish at one hour and 24 hours

(Ivoclar does not contain added calcium)

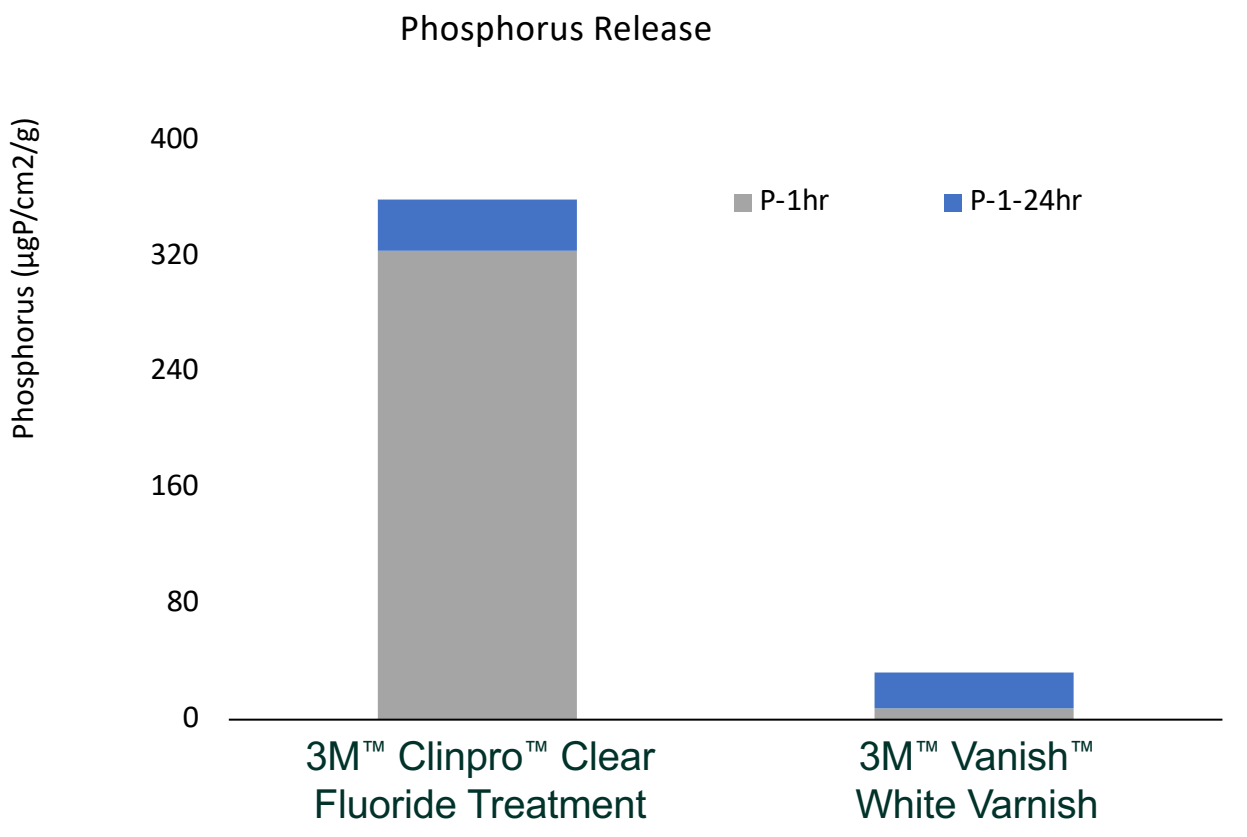- The cumulative phosphorus release of Clinpro Clear Fluoride Treatment is statistically higher than Vanish White Varnish at one hour and 24 hours

(Ivoclar does not contain added phosphorus)

3M™ Vanish™ 5% Sodium Fluoride White Varnish (North America brand) and 3M Clinpro 5% Sodium Fluoride White Varnish (brand outside of North America) are equivalent products.

Source: Ion Release Potential of Different Fluoride-Containing Products. Wang, Y., *J Dent Res* Vol Spec Iss#104:0909